



Installation & Configuration

INDEX

INSTALLATION 3

- quickPX-Installation 3**
- Nikon Camera Control Pro 2 (CCP2)- Installation 4**
- Camera Control Pro 2 - Update 6**

CONFIGURATION..... 7

- Presets..... 7**
- Standard Settings..... 10**
- Camera-specific Adjustments 11**
 - Adjustments for D7000..... 11
 - Adjustments for D7100..... 16
 - Adjustments for D5100..... 21
 - Adjustments for D5300..... 26
 - Adjustments for D300S 31
 - Adjustments for D90..... 35

INSTALLATION

quickPX-Installation


First please make sure, that you are logged in to your pc with an account with administrator rights!

Note: The images show the German version of the windows in the installation process. The English text will differ slightly.

1. Drag the quickPXinstall folder into any partition you wish. Please put it directly into the partition, not in a folder.


Name	Änderungsdatum	Typ	Größe
Buffer	12.01.2015 15:28	Dateiordner	
ProgramData	12.01.2015 15:28	Dateiordner	
Projekte	12.01.2015 15:29	Dateiordner	
quickPX	12.01.2015 15:30	Dateiordner	
Wincmd704	12.01.2015 15:30	Dateiordner	
Windows	12.01.2015 15:30	Dateiordner	
quickPXinstaller	12.01.2015 07:07	Anwendung	13.462 KB

2. Double-click on „quickPXinstaller.exe“.

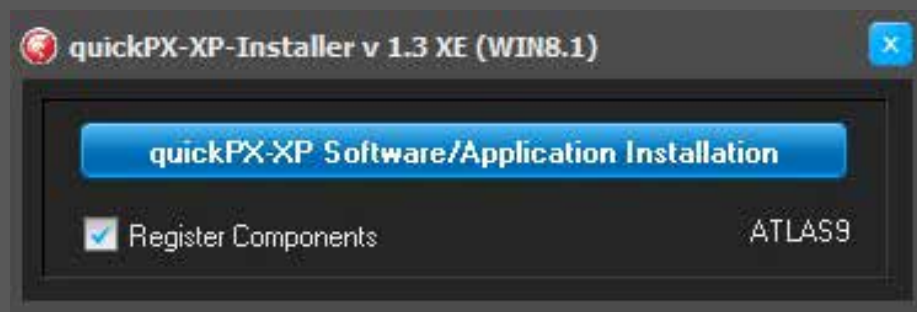


3. Wait until the following notifications are shown and acknowledge them by pressing „OK“.



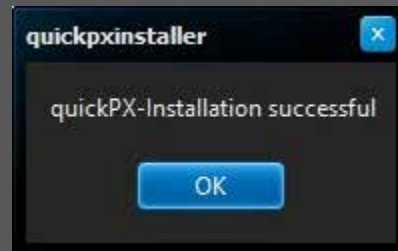

Nikon Camera Control Pro 2 (CCP2)- Installation


1. Open the system partition of your hard drive (in most cases „c“). Then open the folder „program files (x86)“ and the folder „quickPX“.

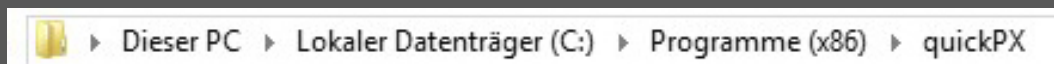

2. Now open the folder „nikon“ and start the program „S-CCPRO_-021300WF-ALLIN-32BIT_.exe“ by double-clicking on it.

Name	Änderungsdatum	Typ	Größe
Nikondeinstall	24.04.2014 13:30	Microsoft Word 9...	16 KB
Regsweeper2	24.04.2014 13:16	Anwendung	15.844 KB
S-CCPRO_-021300WF-ALLIN-32BIT_	31.05.2013 12:22	Anwendung	64.101 KB
S-CCPRO_-021700WU-__EN-32BIT_	24.04.2014 13:26	Anwendung	48.803 KB

3. **IMPORTANT:** Please make sure you choose the language option „English“ in the following window and acknowledge your choice by pressing on the „Next“-button.

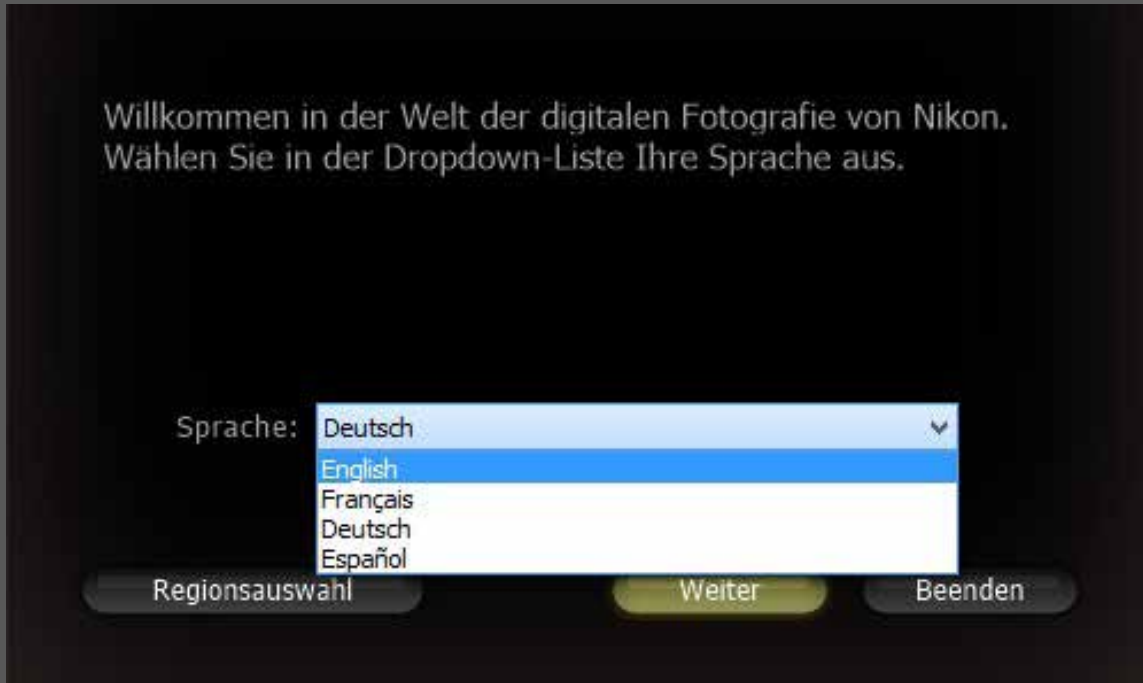

4. Now click on „Install“ in the next window. The InstallShield Wizard will guide you through the installation.



Camera Control Pro 2 - Update

1. Please double-click now on the file „S-CCPRO_-021700WU-___EN-32BIT_“ in the „nikon“ folder off he „quickPXinstall“ folder and acknowledge the update request by clicking on „yes“.



The InstallShield Wizard will now guide you through the installation of the update.

2. As soon as the installation of the update is completed open system control (on Windows 8: right-click on the Windows-symbol in the task bar and then on „system control“) and click on „programs and features“. Search the program „Nikon Message Center 2“ in the list, mark it by left-clicking on it and then click on „Deinstall“ in the command bar above. (Acknowledge your command by pressing the „yes“-button).

Name	Deinstalliert dieses Programm.	Herausgeber	Installiert am	Größe	Version
DriverTurbo		DeskToolsSoft	16.09.2014	12,3 MB	3.2.0
Enfocus PitStop Professional		Enfocus Software	01.10.2014		7.22
gSyncit		Fieldston Software	09.01.2015	13,4 MB	4.0.520
Intel(R) Management Engine Components		Intel Corporation	16.09.2014		10.0.1.1000
Logitech Gaming Software 8.56		Logitech Inc.	02.10.2014	96,6 MB	8.56.109
Microsoft Office 365 Small Business Premium - de-de		Microsoft Corporation	24.12.2014		15.0.4675.1003
Microsoft Office Proofing Tools 2013 - Español		Microsoft Corporation	06.11.2014	22,8 MB	15.0.4420.1017
Microsoft Outlook 2013 - de-de		Microsoft Corporation	24.12.2014		15.0.4675.1003
Microsoft Visual C++ 2008 Redistributable - x64 9.0.3...		Microsoft Corporation	16.09.2014	13,2 MB	9.0.30729
Microsoft Visual C++ 2008 Redistributable - x64 9.0.3...		Microsoft Corporation	01.10.2014	13,2 MB	9.0.30729.6161
Microsoft Visual C++ 2008 Redistributable - x86 9.0.3...		Microsoft Corporation	01.10.2014	10,1 MB	9.0.30729.6161
Microsoft Visual C++ 2010 x64 Redistributable - 10.0...		Microsoft Corporation	07.10.2014	13,8 MB	10.0.40219
Microsoft Visual C++ 2010 x86 Redistributable - 10.0...		Microsoft Corporation	06.10.2014	12,2 MB	10.0.40219
Microsoft Visual C++ 2012 Redistributable (x64) - 11.0...		Microsoft Corporation	07.10.2014	20,5 MB	11.0.61030.0
Microsoft Visual C++ 2012 Redistributable (x86) - 11.0...		Microsoft Corporation	06.10.2014	17,3 MB	11.0.61030.0
Microsoft-Maus- und Tastatur-Center		Microsoft Corporation	19.09.2014		2.3.188.0
Nikon Message Center 2		Nikon	12.01.2015	9,42 MB	2.1.0
Norton Internet Security		Symantec Corporation	16.09.2014		21.6.0.32
NVIDIA 3D Vision Treiber 337.88		NVIDIA Corporation	16.09.2014		337.88
NVIDIA Grafiktreiber 337.88		NVIDIA Corporation	16.09.2014		337.88

Nikon Produktversion: 2.1.0

CONFIGURATION

Presets

1. Please turn on the camera first by turning the switch on the right at the top of the camera to „On“.
2. Now open the program „Nikon Camera Control Pro 2“ (CCP2). (Note: When starting the program for the first time you'll be requested to type in a licence key. You can find this key on the CCP2 box. With a click on „Trial“ you can begin a trial period of 30 days. During this phase you can use the program without restrictions. In this case please make sure that you always start up the CCP2 program before quickPX to avoid malfunctions.)
3. Now click on the item „Tools“ on the menu bar. A pop-up menu opens itself. Choose the menu item „Transfer Options“ with a left-click.
4. In the following window there are several adjustments to make:

Transfer Options

Destination folder:
C:\Users\Public\quickpx\CamInFolder Browse...

File name to be used:
Img002133 Edit...
This setting will not apply if images are recorded to both the memory card and PC.

When a new image is received from the camera:
Do nothing


Attach IPTC Information
Note (D4 and D4S only): If the IPTC data was embedded in the camera, the IPTC information that is set on this dialog will not be attached.
IPTC Information...

Copy Shooting Data to IPTC Caption

Embed ICC Profile


OK Cancel Help

- a. Set the destination folder as depicted in the picture above. (system-partition:\Users\Public\quick-px\CamInFolder)
- b. Click on the button „Edit...“ next to the menu item „File name to be used:“. In the following window make sure that the length of the number in the file names are set on „6 digits“. Afterwards click on „OK“.

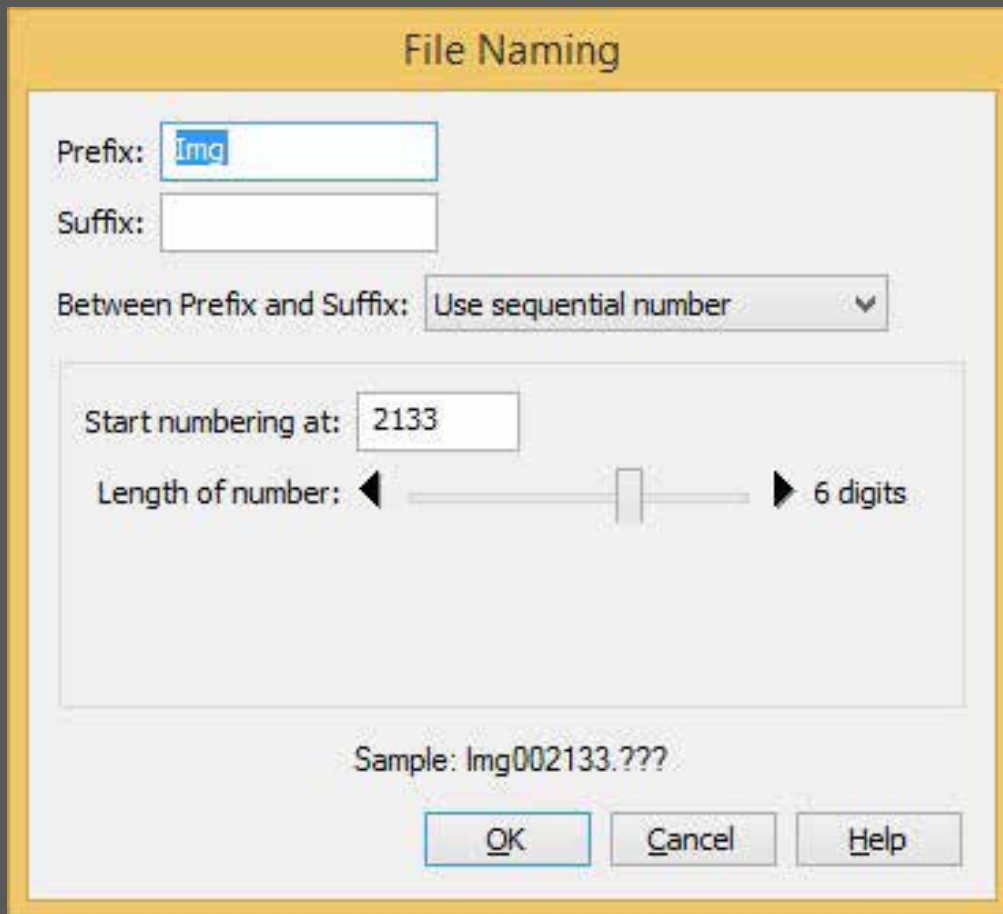


- c. Make sure, that the boxes next to the phrases „Attach IPTC Information“, „Copy Shooting Data to IPTC Caption“ and „Embed ICC Profile“ are all checked
- d. To conclude this step click on „OK“.
5. Now once again click on „Tools“ on the menu bar and then on „Options“ in the following pop-up menu.

6. In the following window click on the card „Color Management“. Click on the button „Browse...“.

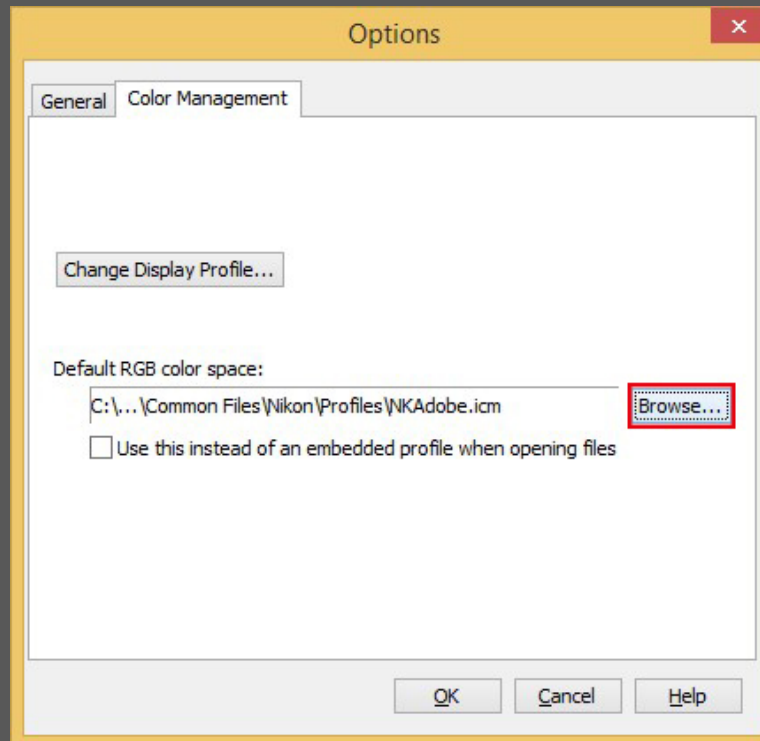

7. Now mark the profile „NKAdobe“ with a left click and click on „Open“.

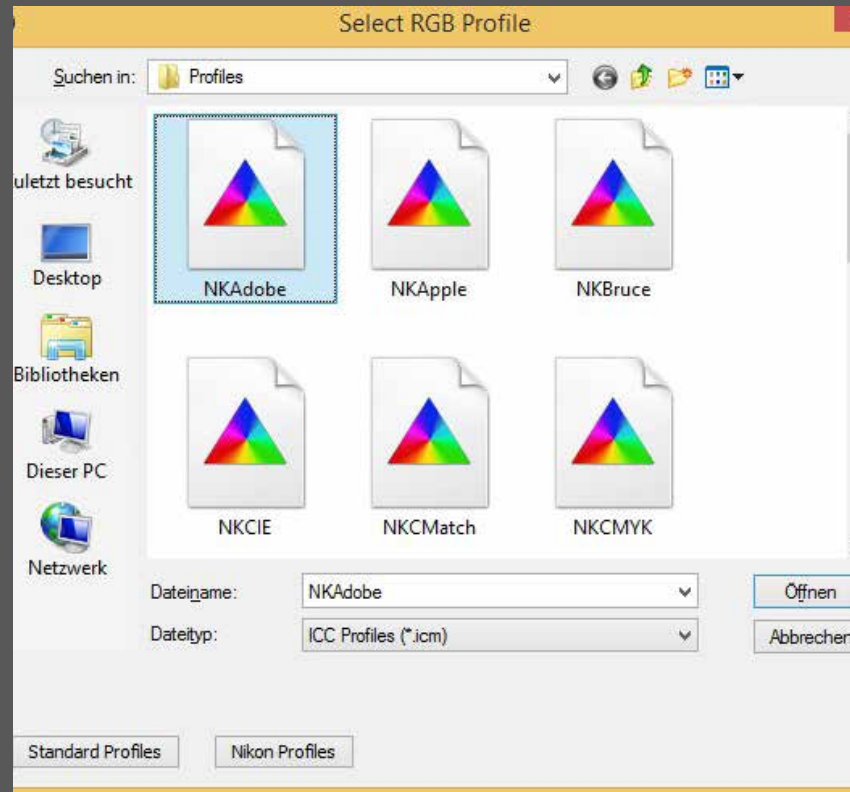


8. Now close the „Options“-window by pressing „OK“.


Standard Settings

NOTE: The menus and interfaces may differ from those depicted inside this manual depending on the software version of CCP2, the type of the used camera and it's intern firmware version.

1. Please click on the button „Set“ in the top are on the left of the quickPX interface.

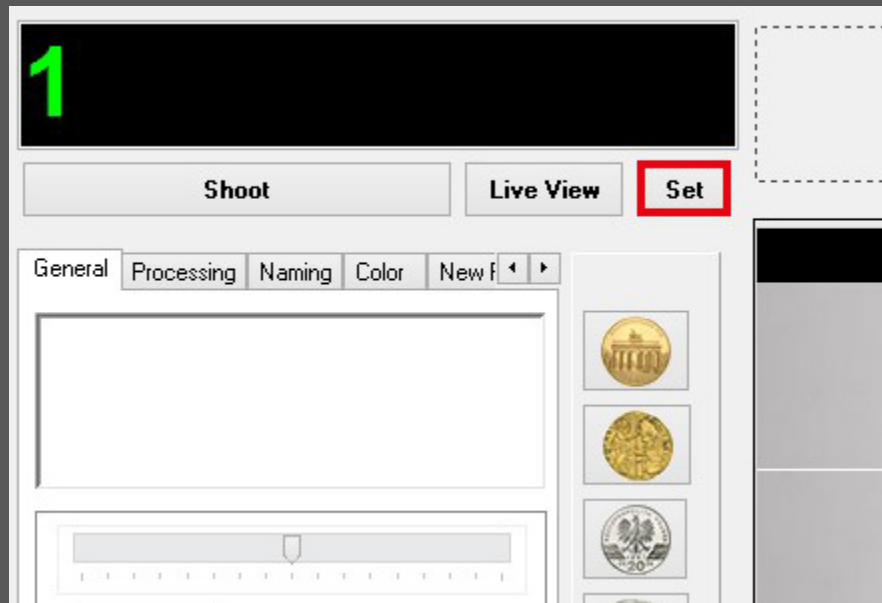

2. The feature „Enable Controls on Camera Body“ is activated by default. To enable quickPX to take over the control off he camera it must be deactivated. Fort hat purpose click on the item „Camera“ in the tast bar of Camera Control Pro 2 and then in the context menu on the item „Enable Controls on Camera Body“, so the check on the left of it disappears.

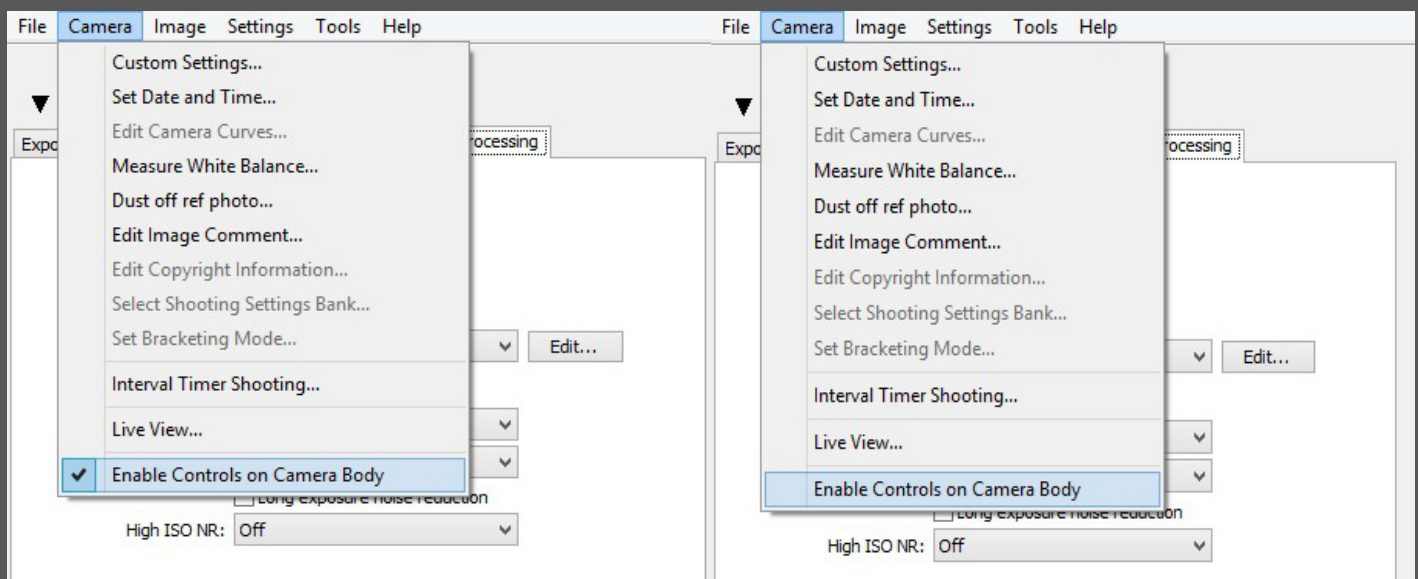


Camera-specific Adjustments

Adjustments for D7000


NOTE: It is necessary, that you adjust time and date on your camera manually when bringing it to service for the first time. Otherwise the CCP2 will not detect the camera!

NOTE: Please ensure that the rotary switch, located on the topside of the camera on the left, is set on „M“ (Manual), so quickPX can take control over the camera.

1. Please open quickPX at first. (CCP2 will also be started automatically. NOTE: If you are using the trial version of CCP2, please start CCP2 before quickPX!)
2. First please click on the card „Presets“ in the adjustments area on the left of the quickPX interface. You may have to use the arrow-buttons on the right of the card bar to unveil it. Now ensure that the item „Switch live view before shooting“ is unchecked.

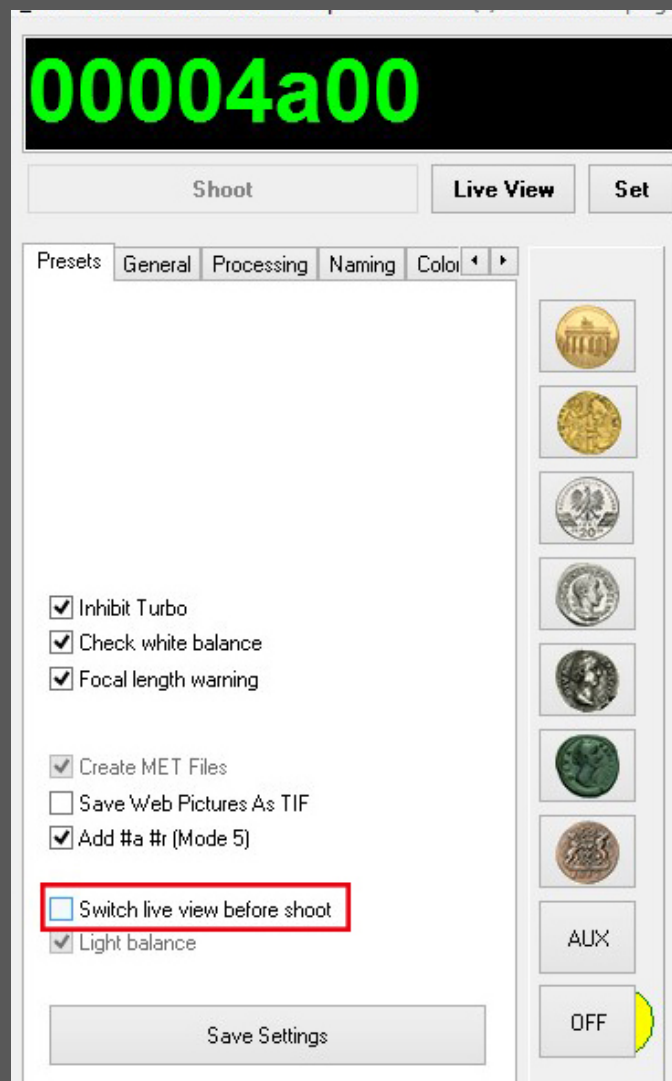


3. To make the remaining adjustments please click on the button „Set“ in the top area of the interface.

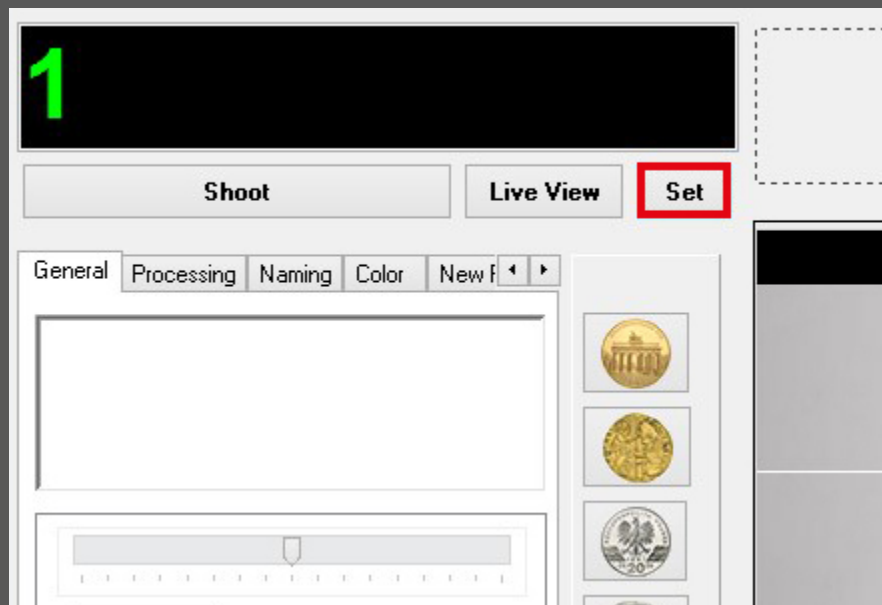


4. Now please ensure that all settings in categories in this menu match the settings depicted in the following pictures.


Exposure 1:

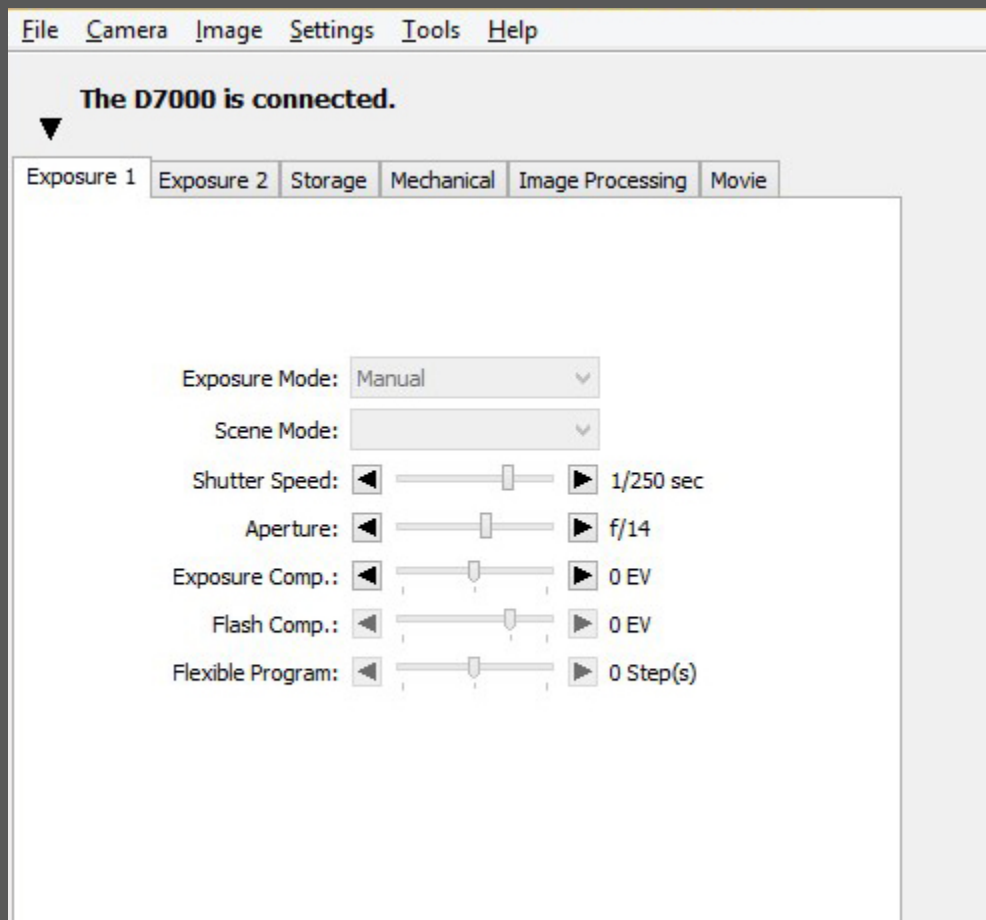

Exposure 2:

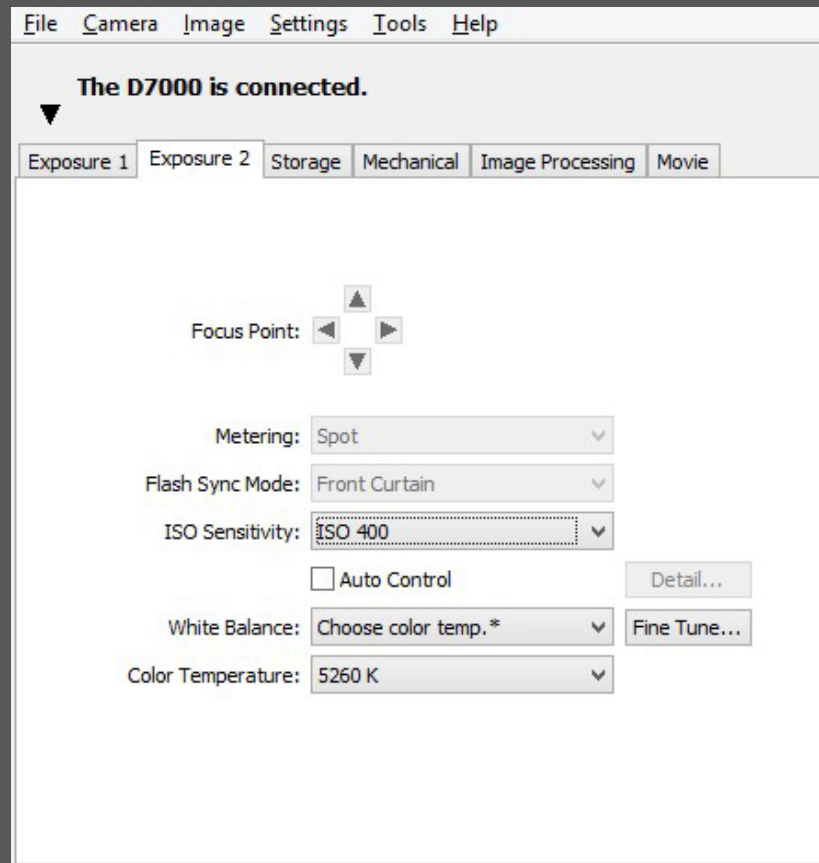

Storage:

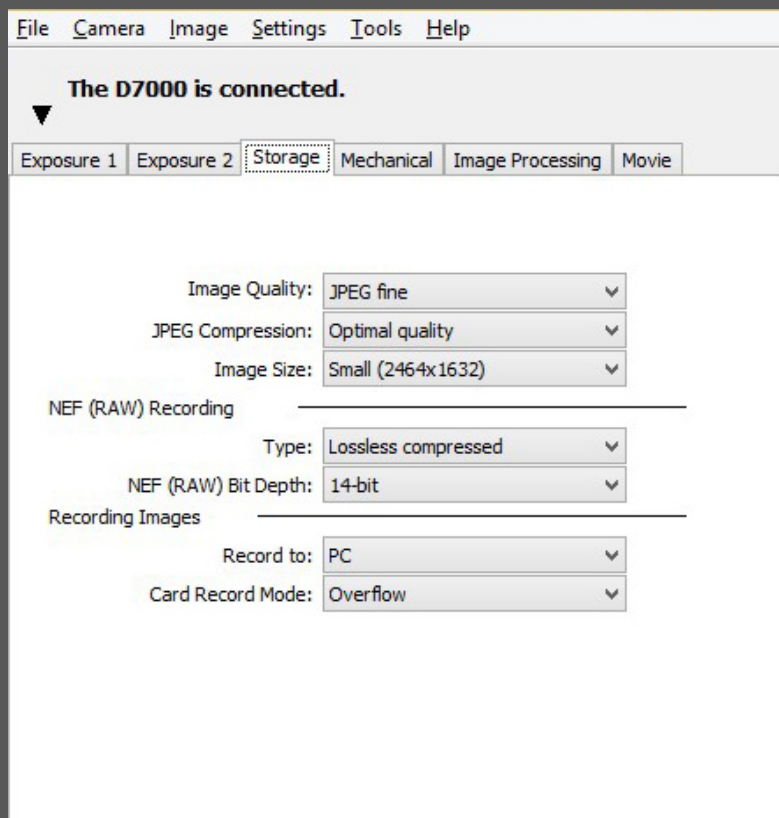

Mechanical:



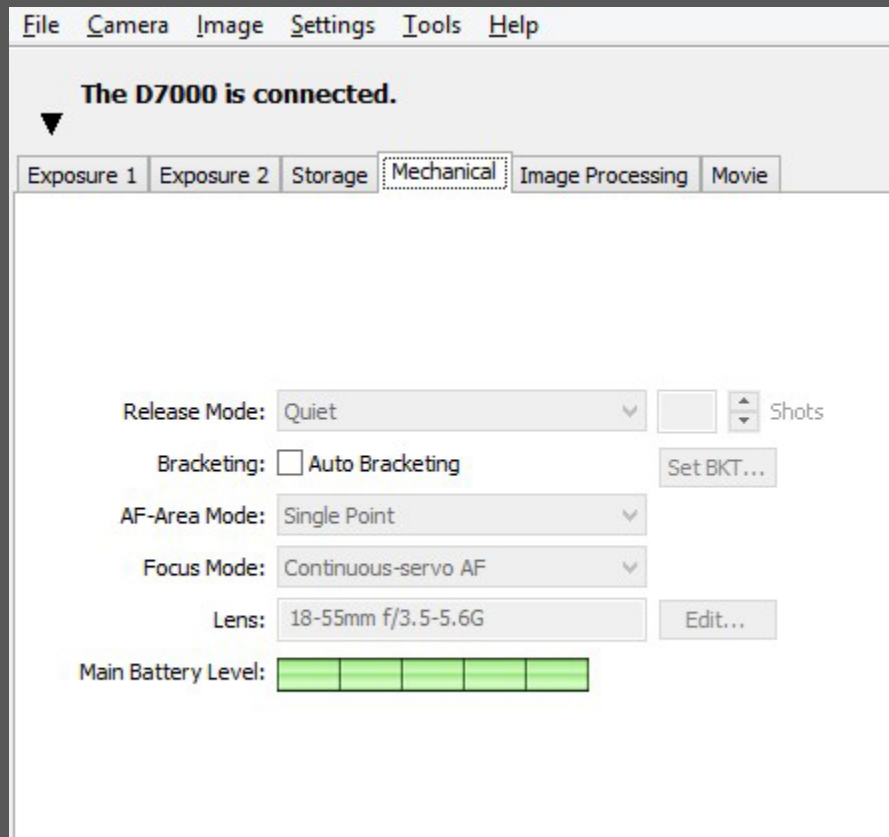


Image Processing:

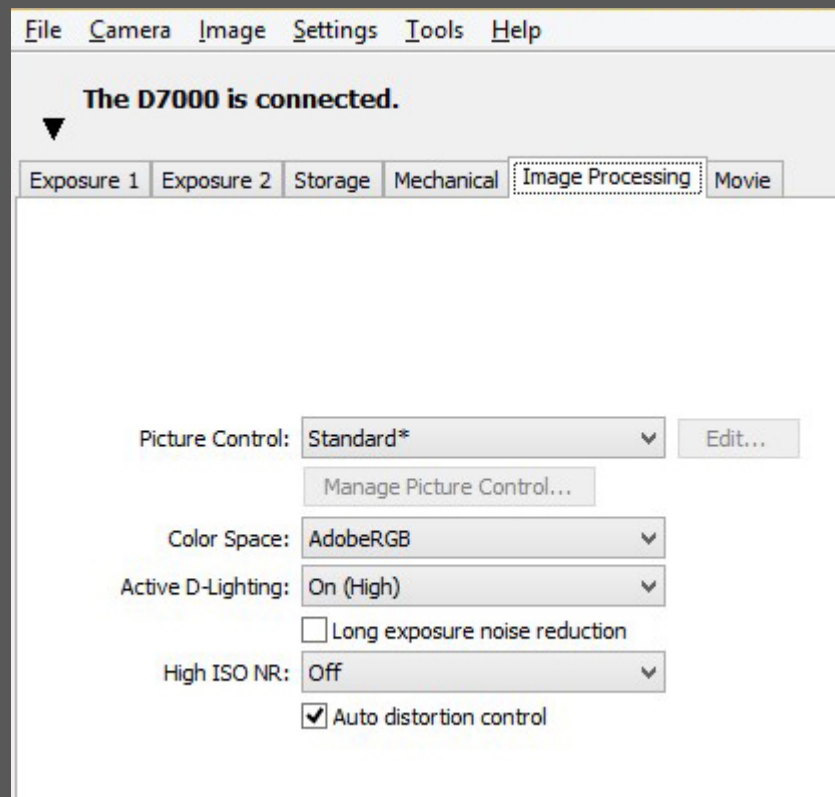

Movie:

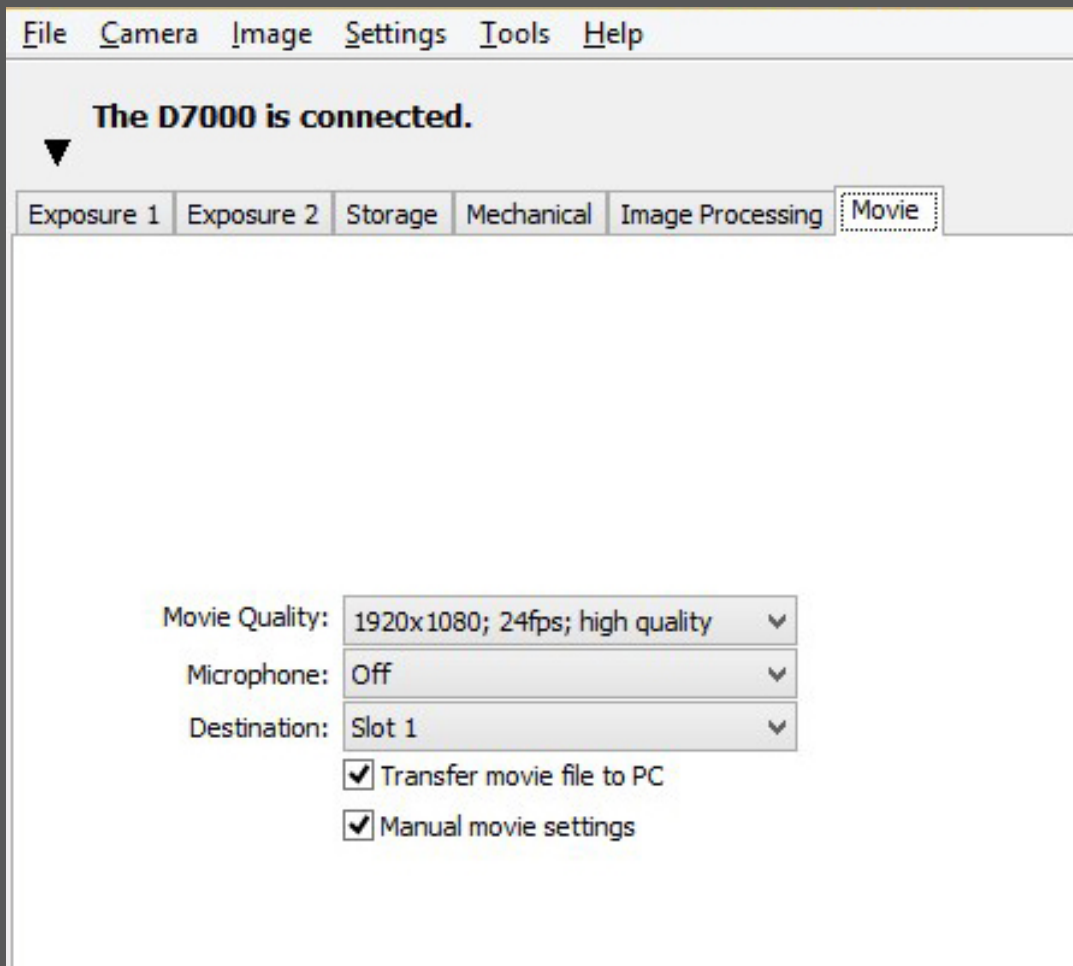


Adjustments for D7100


NOTE: It is necessary, that you adjust time and date on your camera manually when bringing it to service for the first time. Otherwise the CCP2 will not detect the camera!

NOTE: Please ensure that the rotary switch, located on the topside of the camera on the left, is set on „M“ (Manual), so quickPX can take control over the camera.

1. Please open quickPX at first. (CCP2 will also be started automatically. NOTE: If you are using the trial version of CCP2, please start CCP2 before quickPX!)
2. First please click on the card „Presets“ in the adjustments area on the left of the quickPX interface. You may have to use the arrow-buttons on the right of the card bar to unveil it. Now ensure that the item „Switch live view before shooting“ is unchecked.

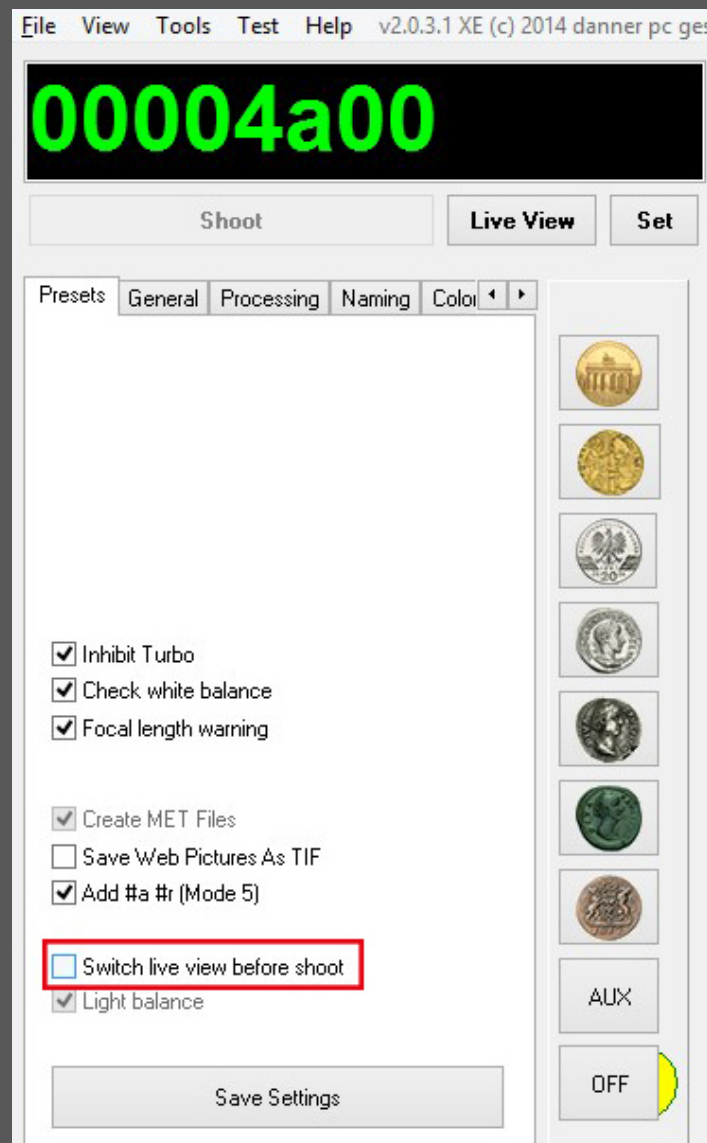


- To make the remaining adjustments please click on the button „Set“ in the top area of the interface.

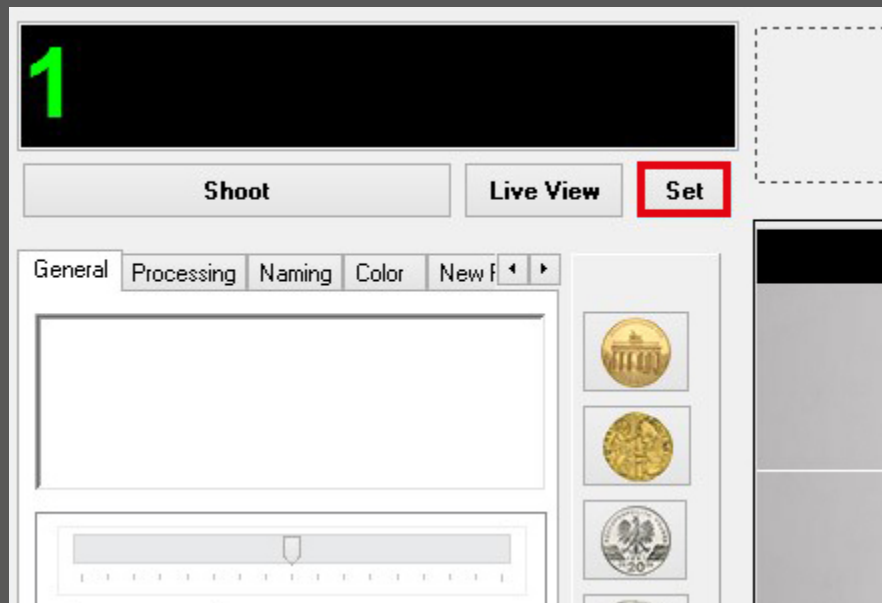


- Now please ensure that all settings in categories in this menu match the settings depicted in the following pictures.


Exposure 1:

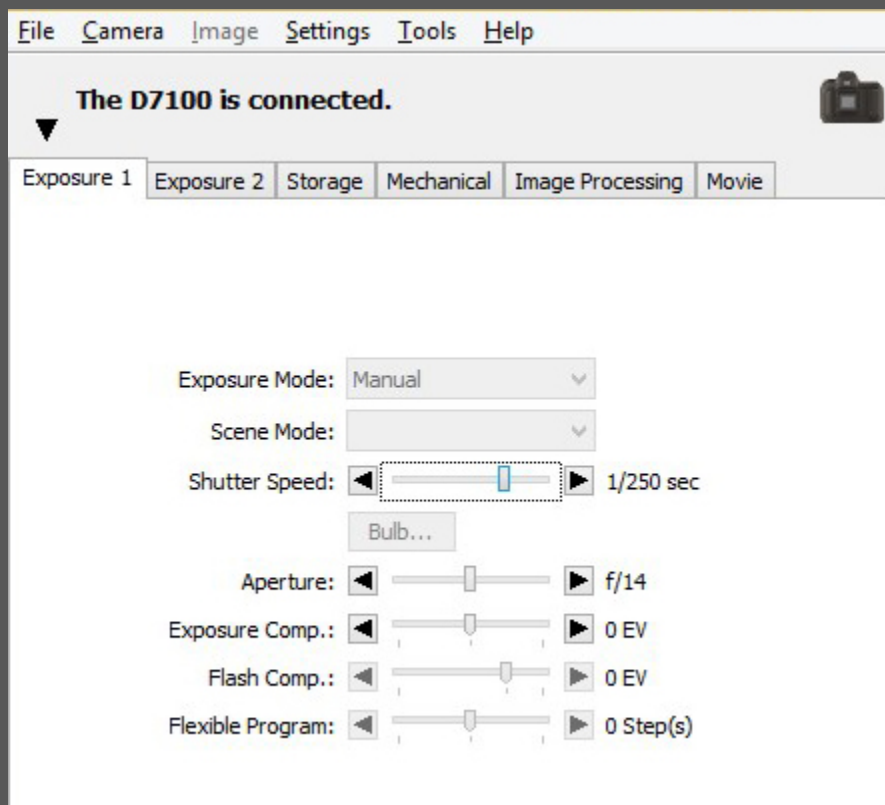

Exposure 2:

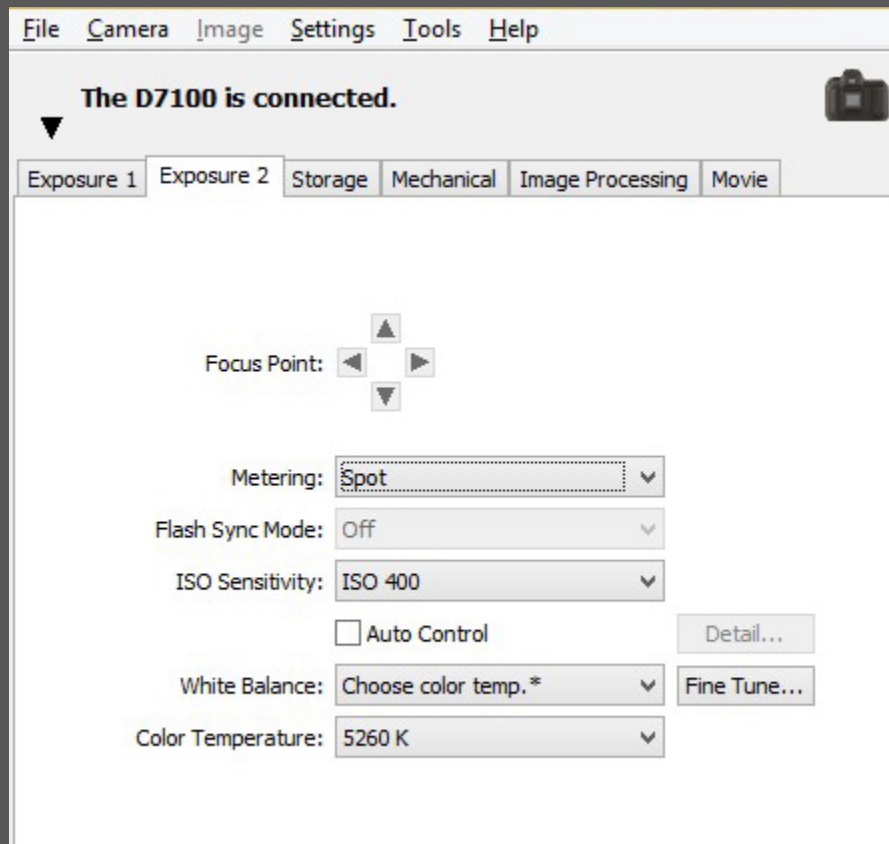

Storage:

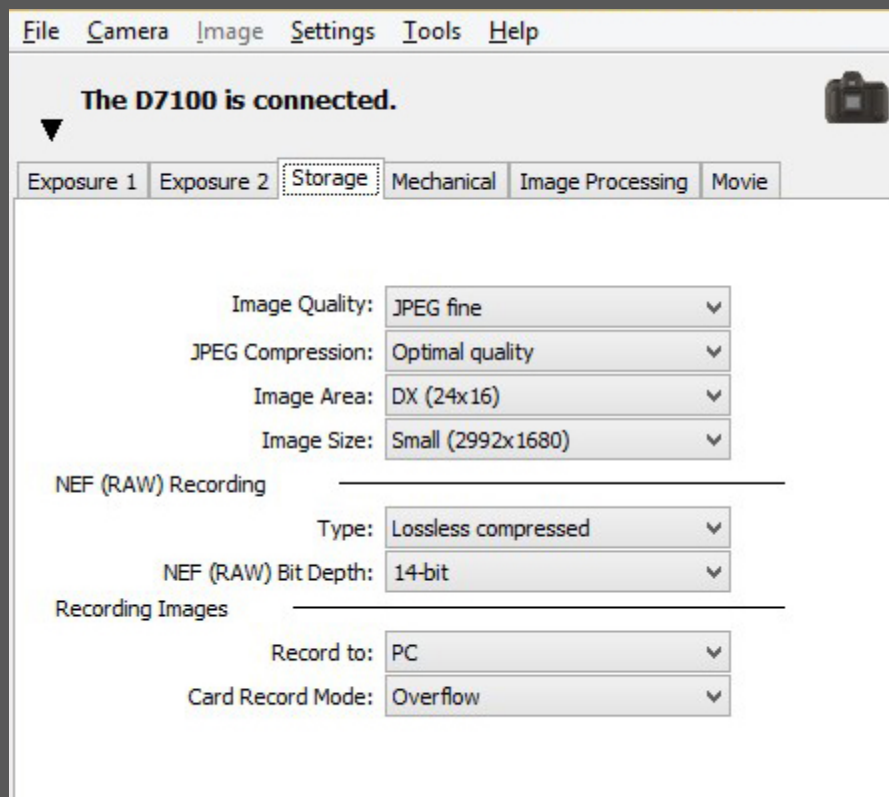

Mechanical:



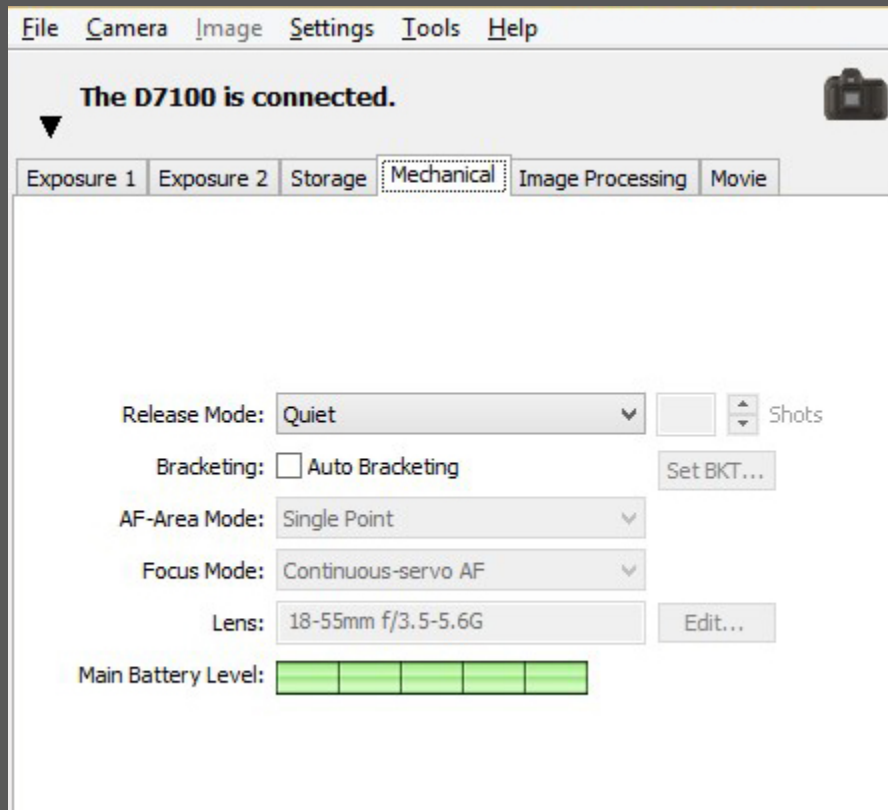


Image Processing:

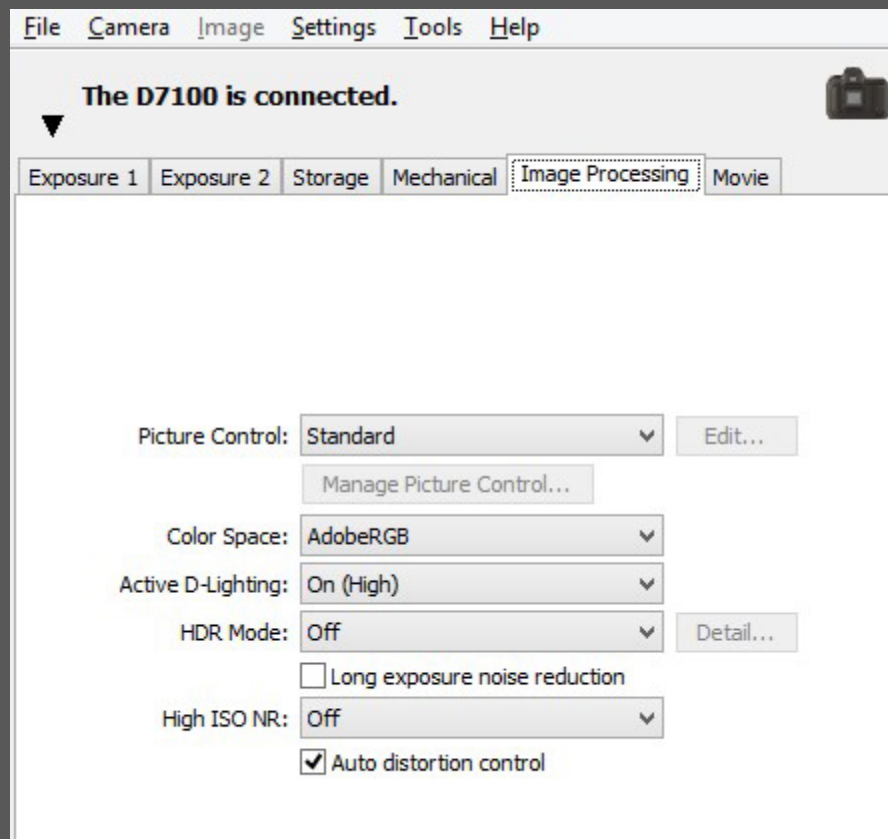

Movie:

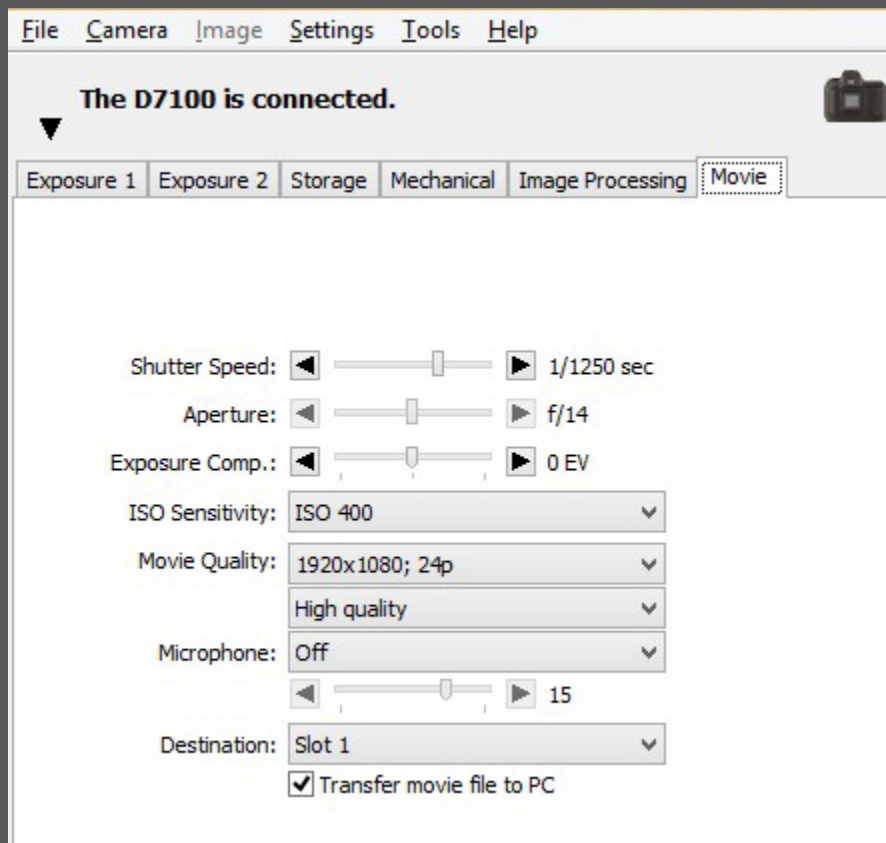


Adjustments for D5100


NOTE: For the D5100 it is necessary, that you open the hinged display and adjust time and date on your camera manually when bringing it to service for the first time. Otherwise the CCP2 will not detect the camera!

NOTE: Please ensure that the rotary switch, located on the topside of the camera on the right, is set on „M“ (Manual), so quickPX can take control over the camera.

1. Please open quickPX at first. (CCP2 will also be started automatically. **NOTE:** If you are using the trial version of CCP2, please start CCP2 before quickPX!)
2. First please click on the card „Presets“ in the adjustments area on the left of the quickPX interface. You may have to use the arrow-buttons on the right of the card bar to unveil it. Now ensure that the item „Switch live view before shooting“ is unchecked.

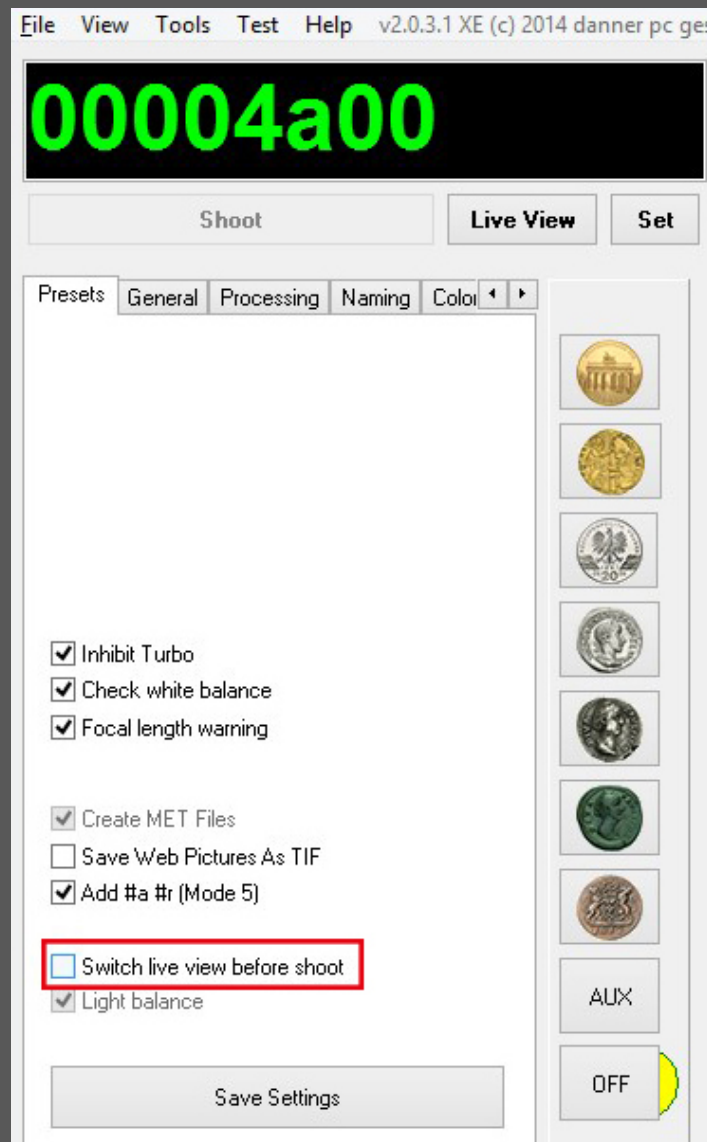


3. To make the remaining adjustments please click on the button „Set“ in the top area of the interface.

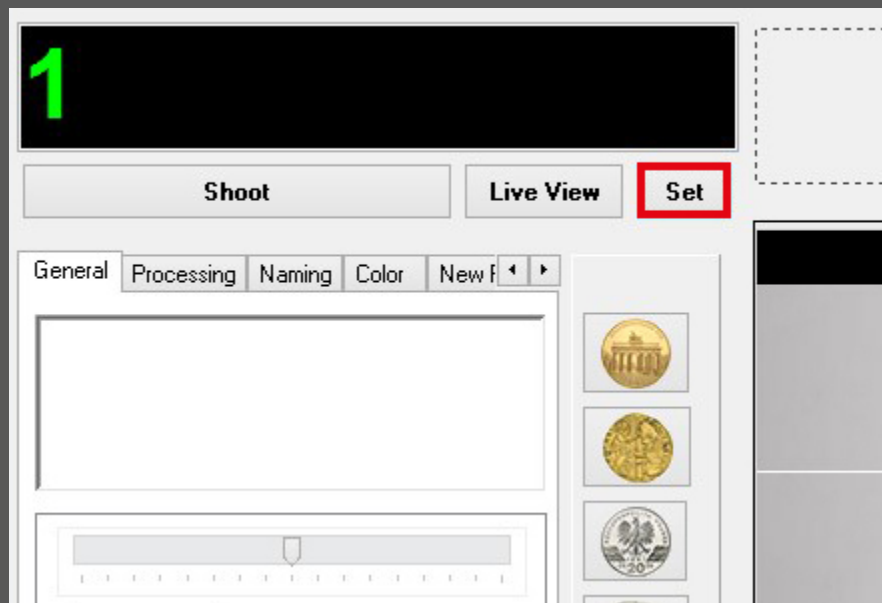


4. Now please ensure that all settings in categories in this menu match the settings depicted in the following pictures.


Exposure 1:

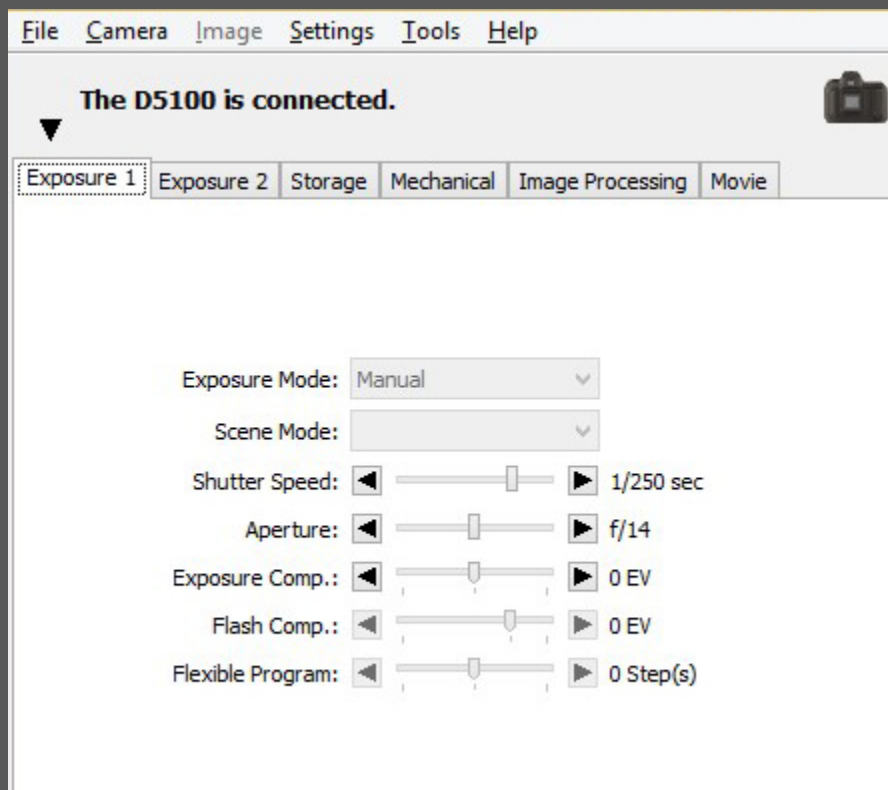

Exposure 2:

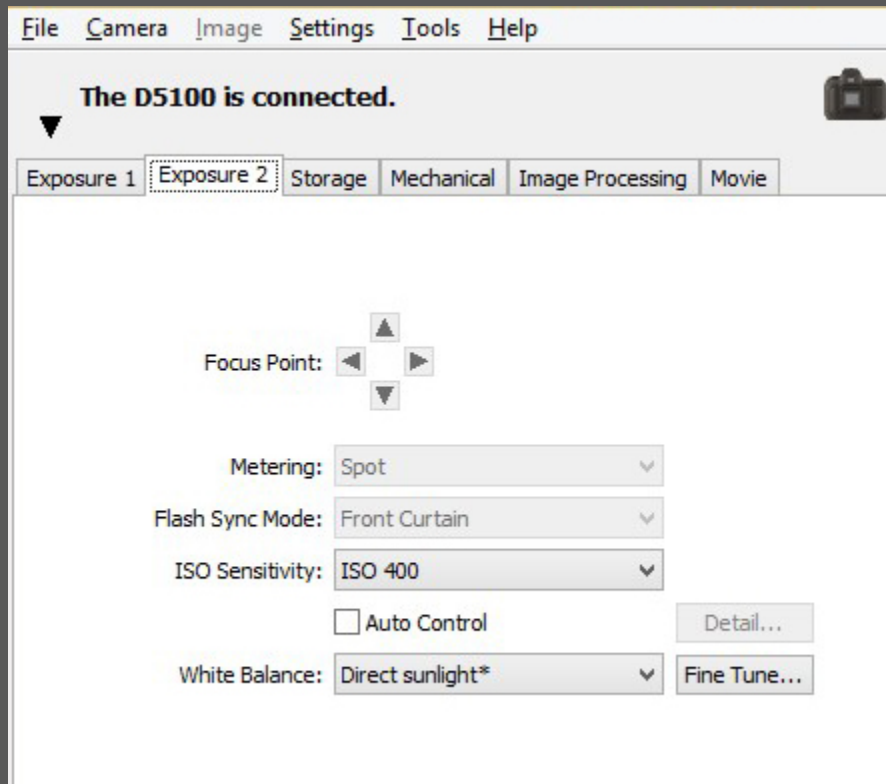

Storage:

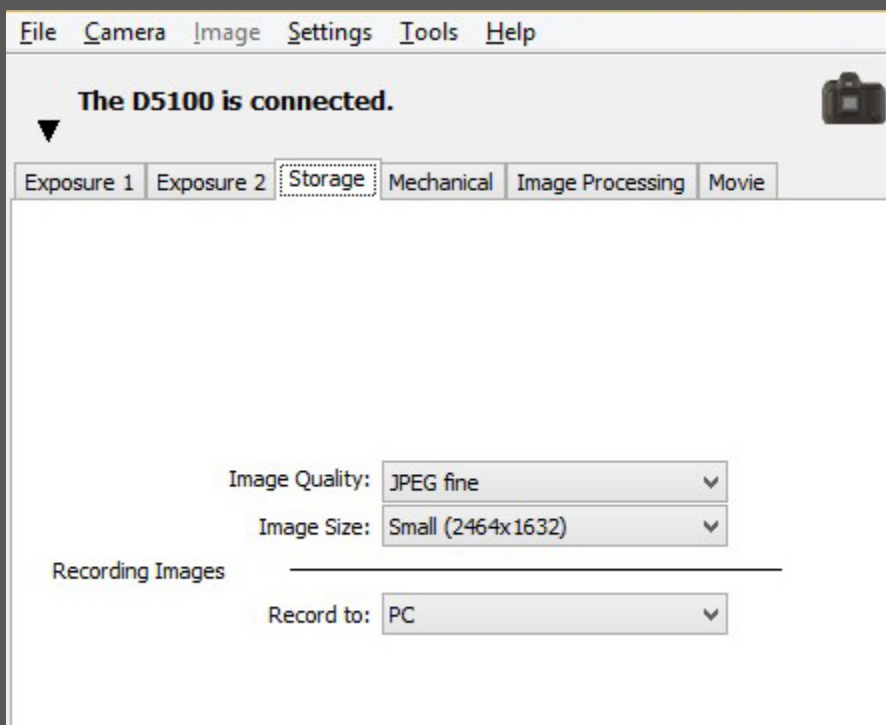

Mechanical:



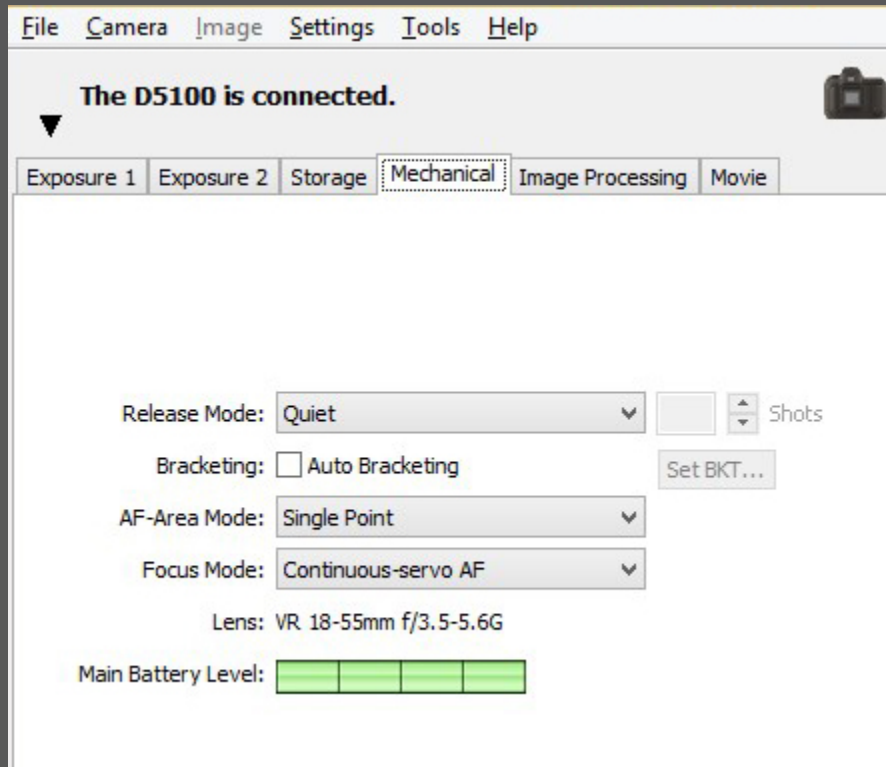


Image Processing:

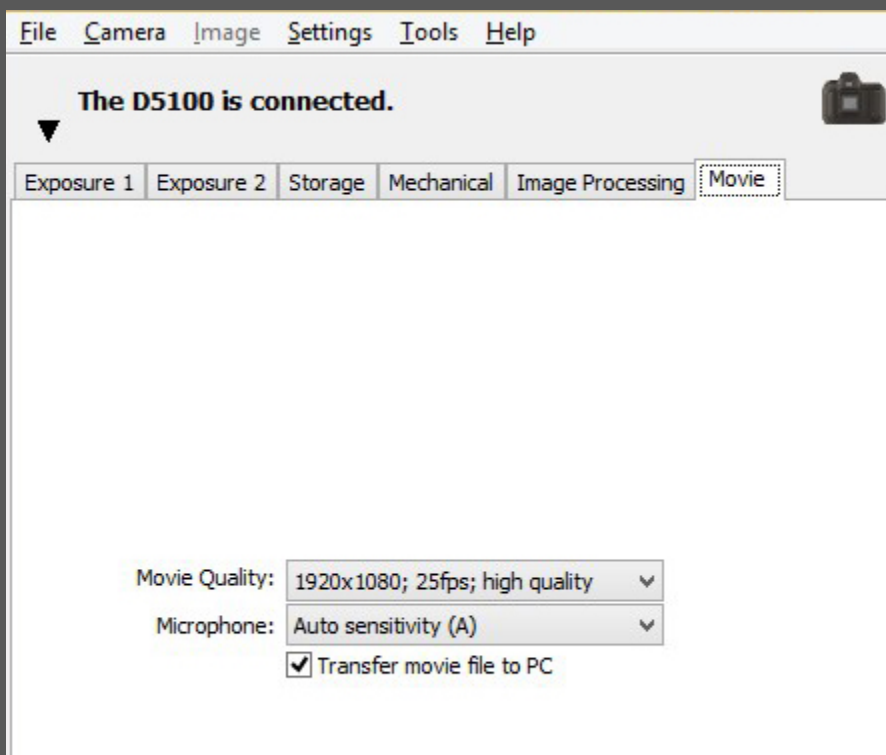

Movie:

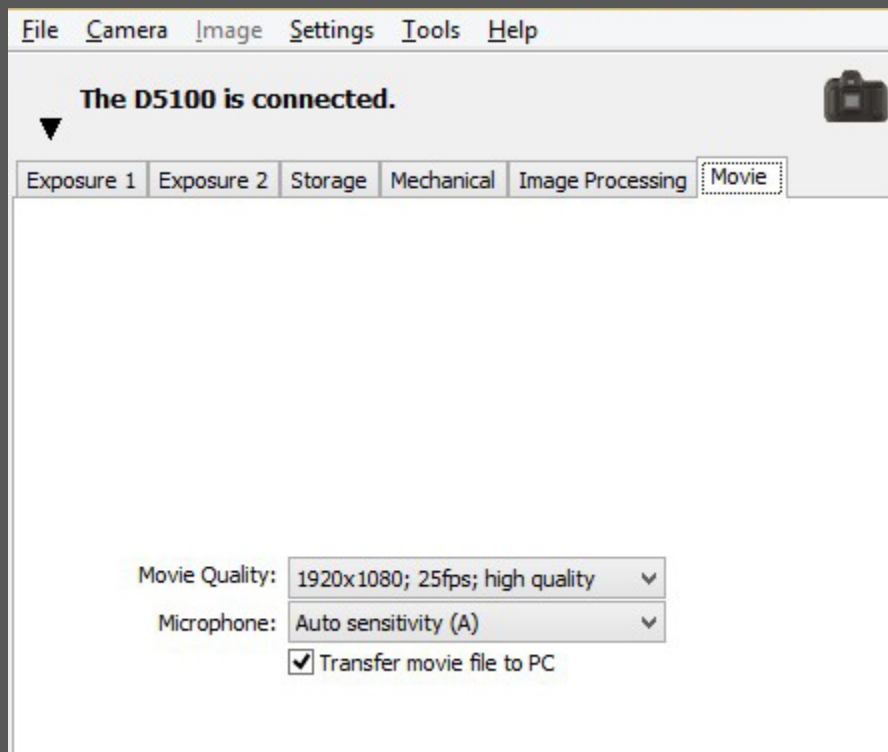


Adjustments for D5300


NOTE: It is necessary, that you adjust time and date on your camera manually when bringing it to service for the first time. Otherwise the CCP2 will not detect the camera!

NOTE: Please ensure that the rotary switch, located on the topside of the camera on the right, is set on „M“ (Manual), so quickPX can take control over the camera.

1. Please open quickPX at first. (CCP2 will also be started automatically. NOTE: If you are using the trial version of CCP2, please start CCP2 before quickPX!)
2. First please click on the card „Presets“ in the adjustments area on the left of the quickPX interface. You may have to use the arrow-buttons on the right of the card bar to unveil it. Now ensure that the item „Switch live view before shooting“ is unchecked.

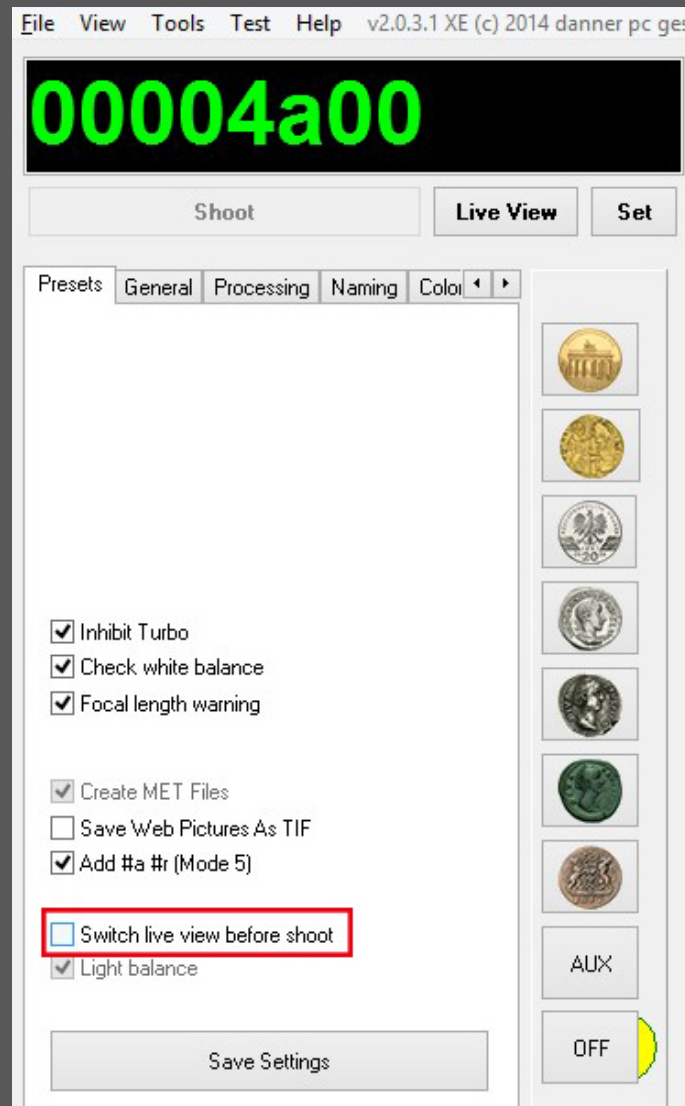


3. To make the remaining adjustments please click on the button „Set“ in the top area of the interface.

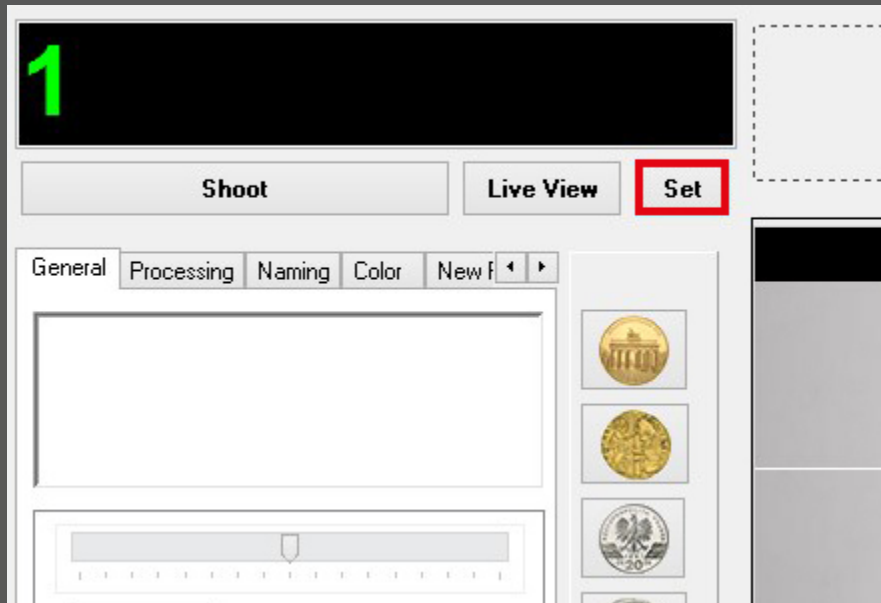


4. Now please ensure that all settings in categories in this menu match the settings depicted in the following pictures.


Exposure 1:

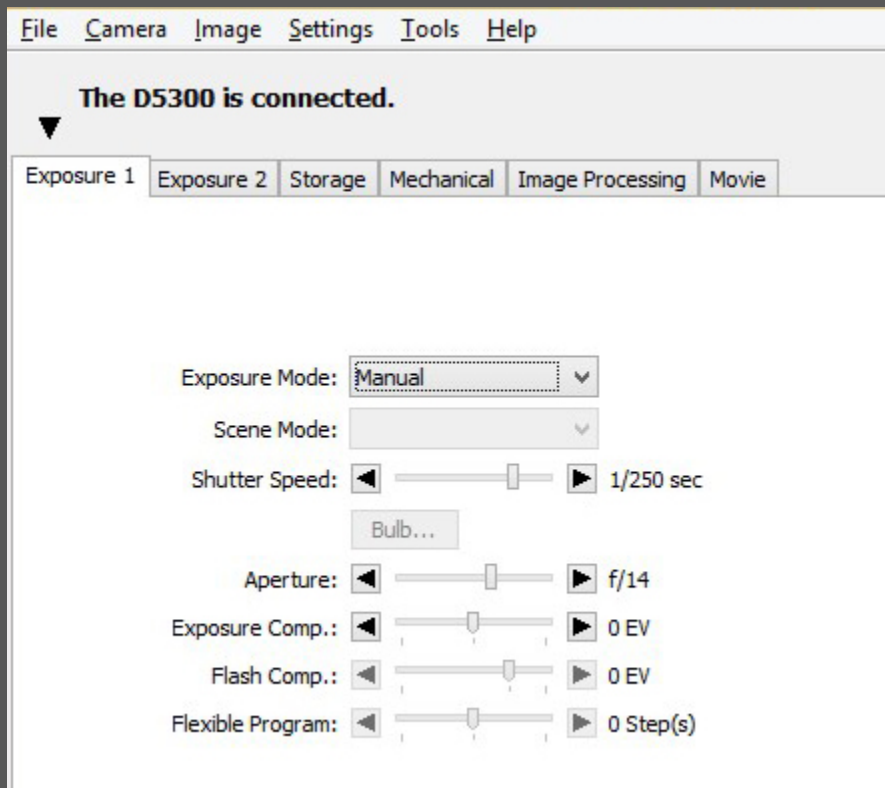

Exposure 2:

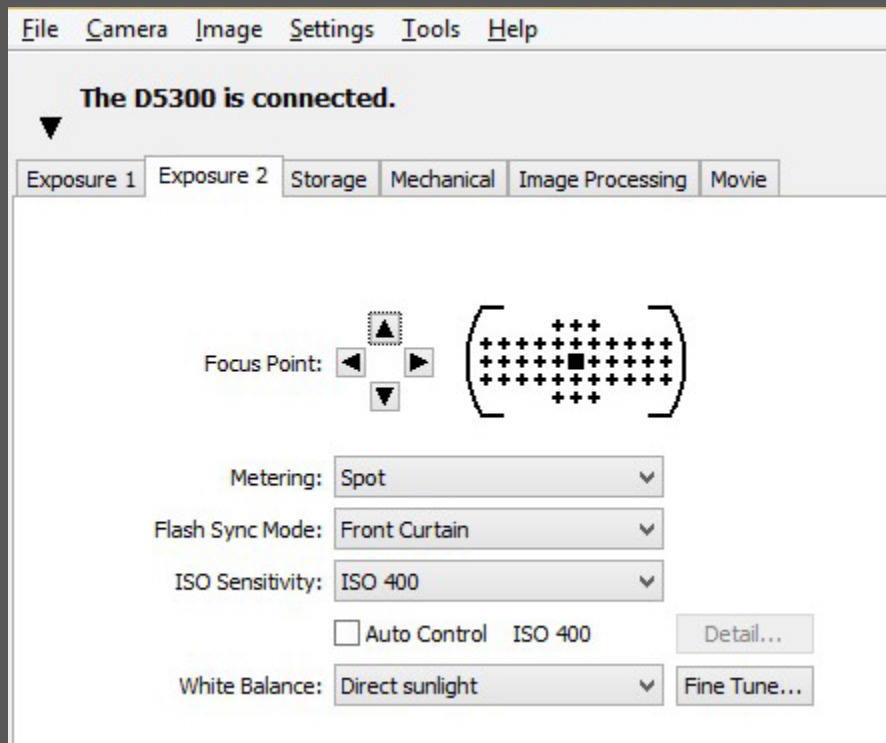

Storage:

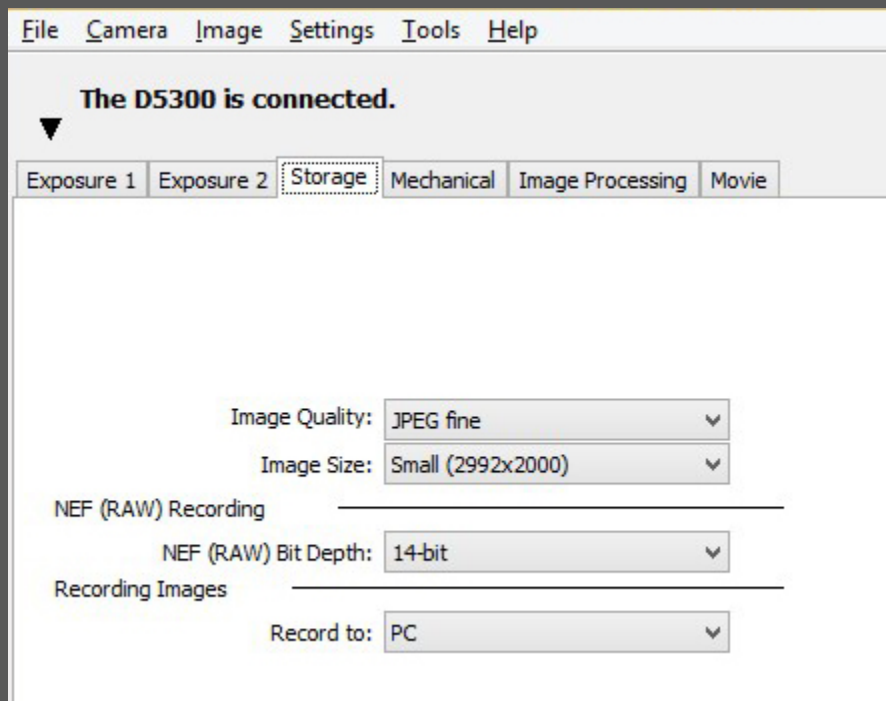

Mechanical:



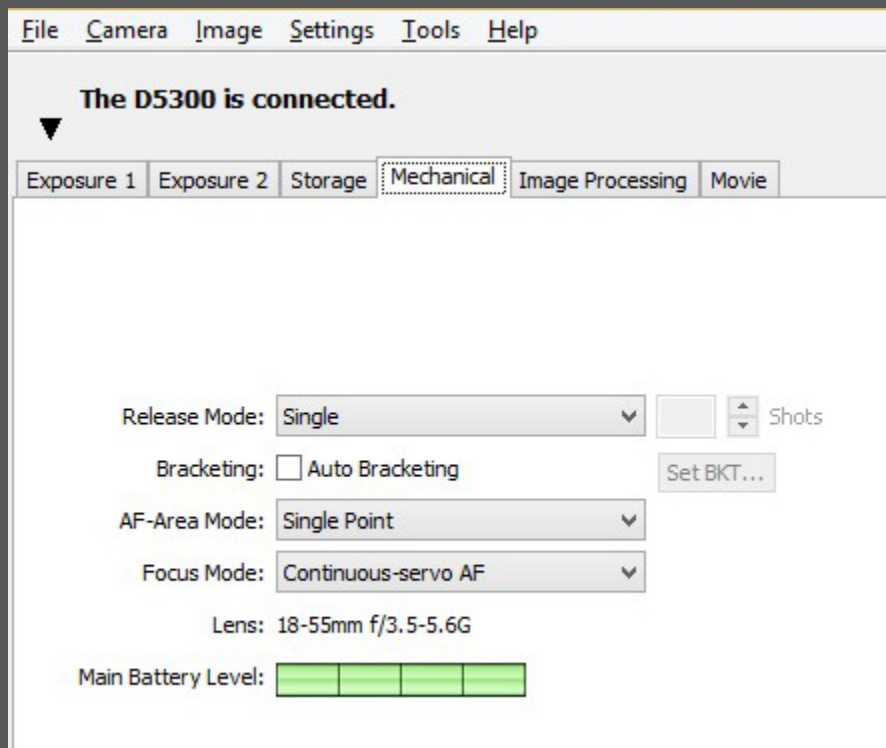


Image Processing:

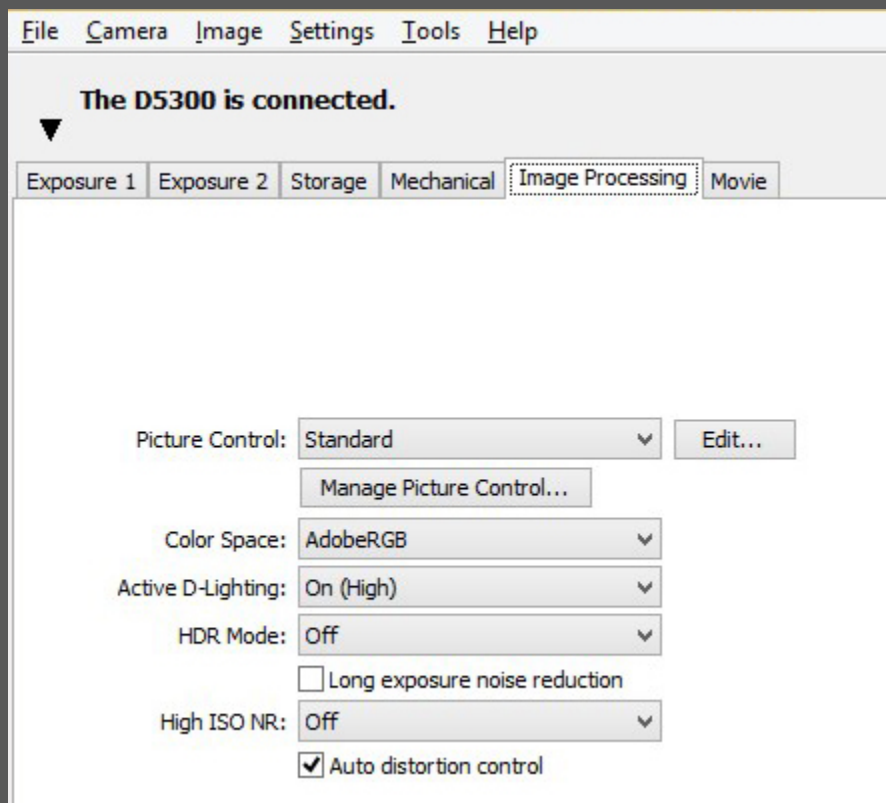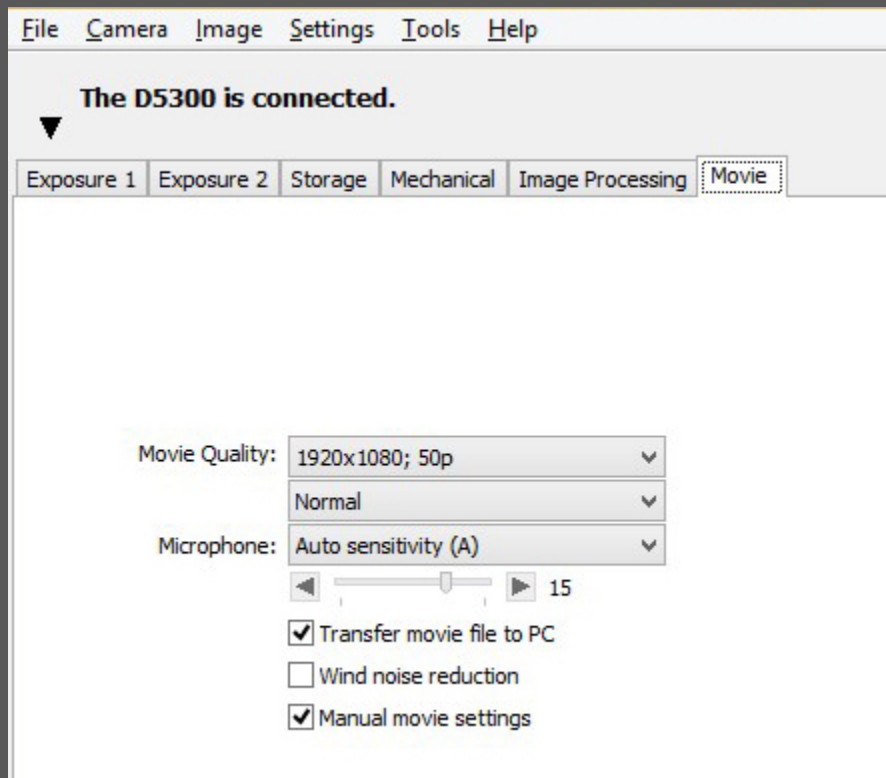


Movie:

Adjustments for D300S


NOTE: It is necessary, that you adjust time and date on your camera manually when bringing it to service for the first time. Otherwise the CCP2 will not detect the camera!

NOTE: Please ensure that you set the photography mode on your camera on „M“ (Manual), so quickPX can take over control over the camera.

1. Please open quickPX at first. (CCP2 will also be started automatically. NOTE: If you are using the trial version of CCP2, please start CCP2 before quickPX!)
2. First please click on the card „Presets“ in the adjustments area on the left of the quickPX interface. You may have to use the arrow-buttons on the right of the card bar to unveil it. Now ensure that the item „Switch live view before shooting“ is checked.

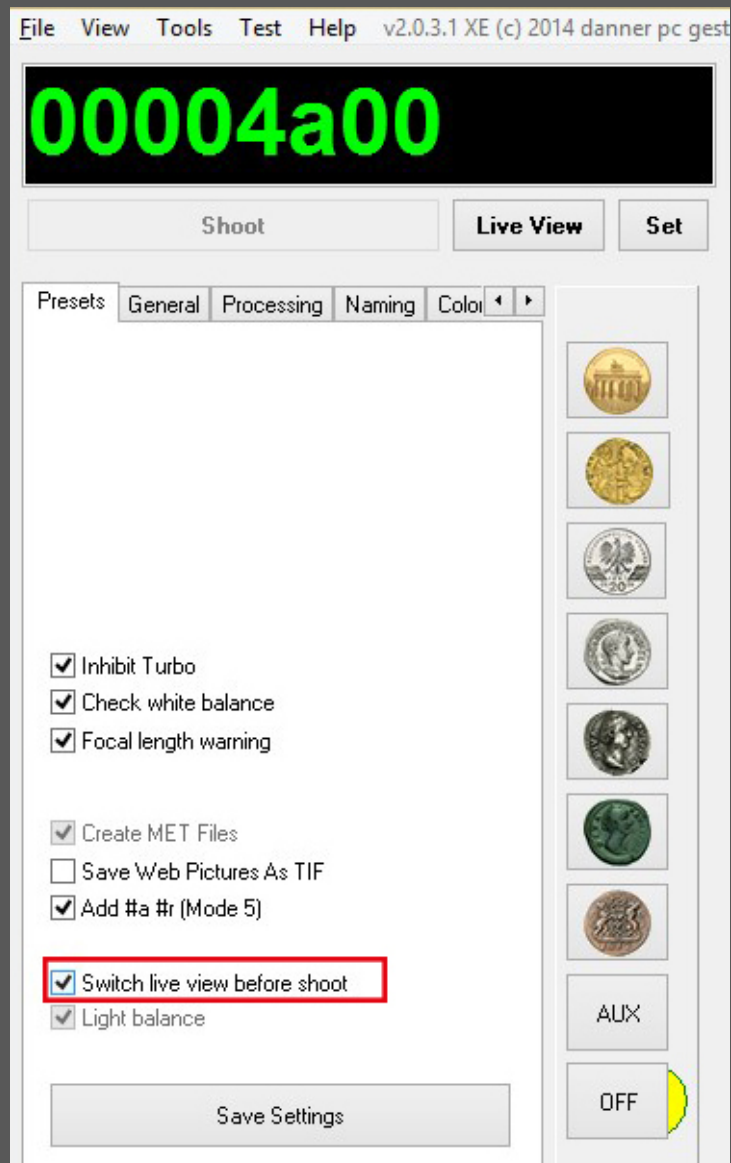


3. To make the remaining adjustments please click on the button „Set“ in the top area of the interface.

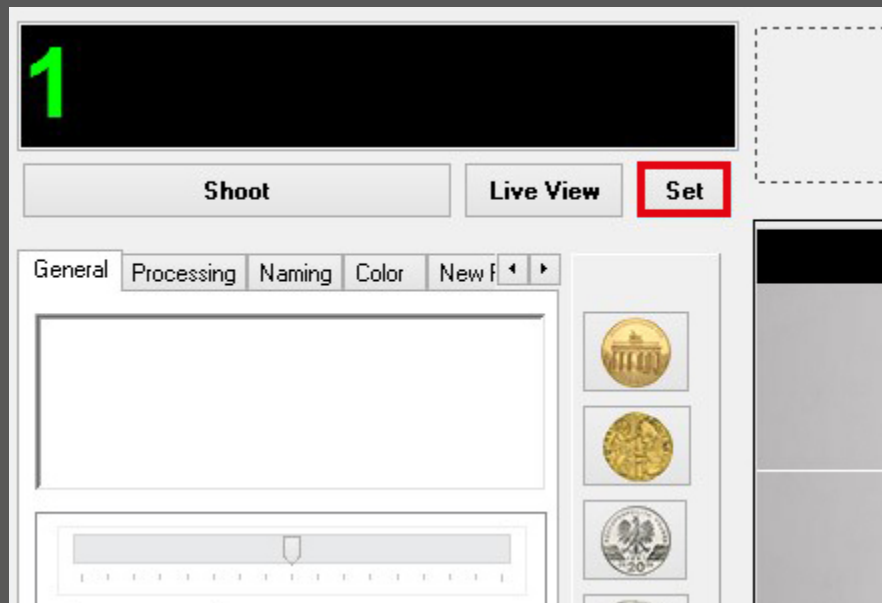


4. Now please ensure that all settings in categories in this menu match the settings depicted in the following pictures.


Exposure 1:

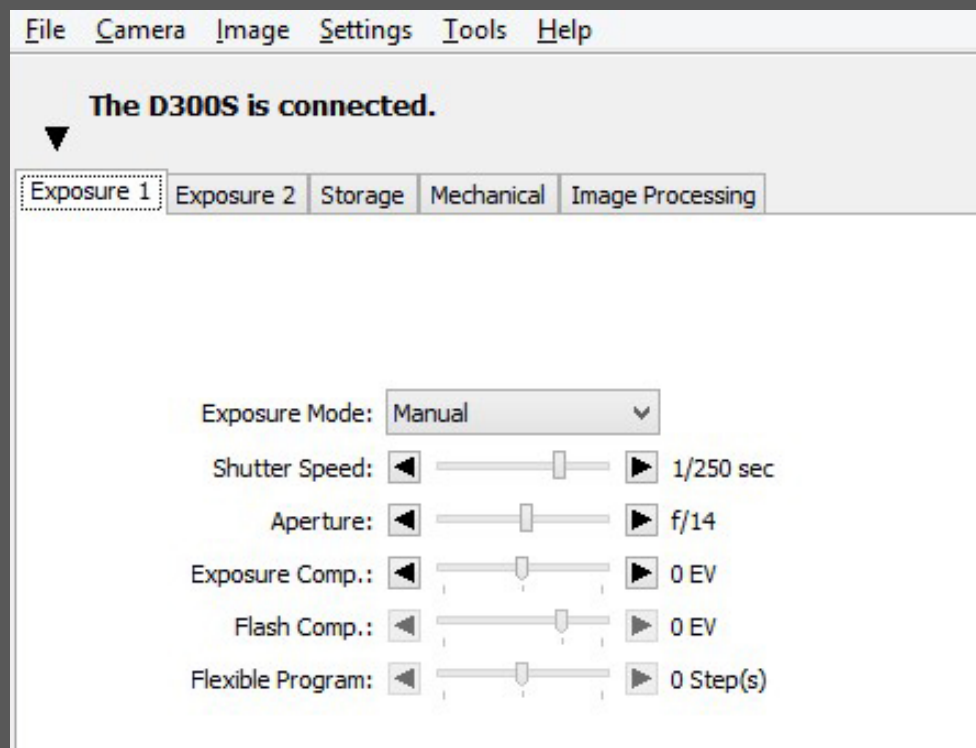

Exposure 2:

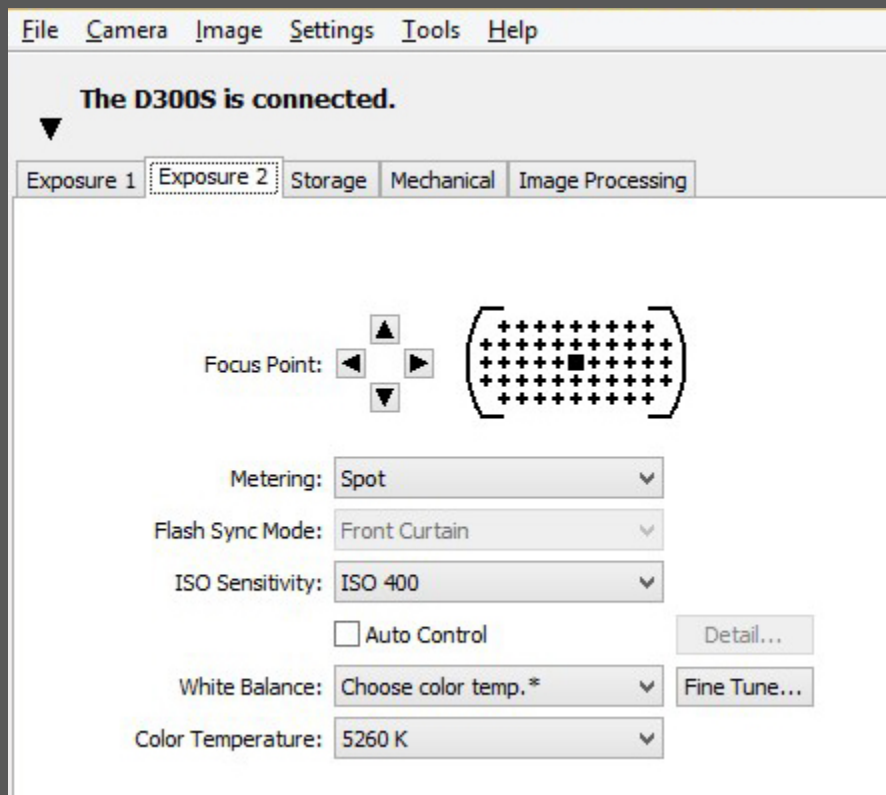

Storage:

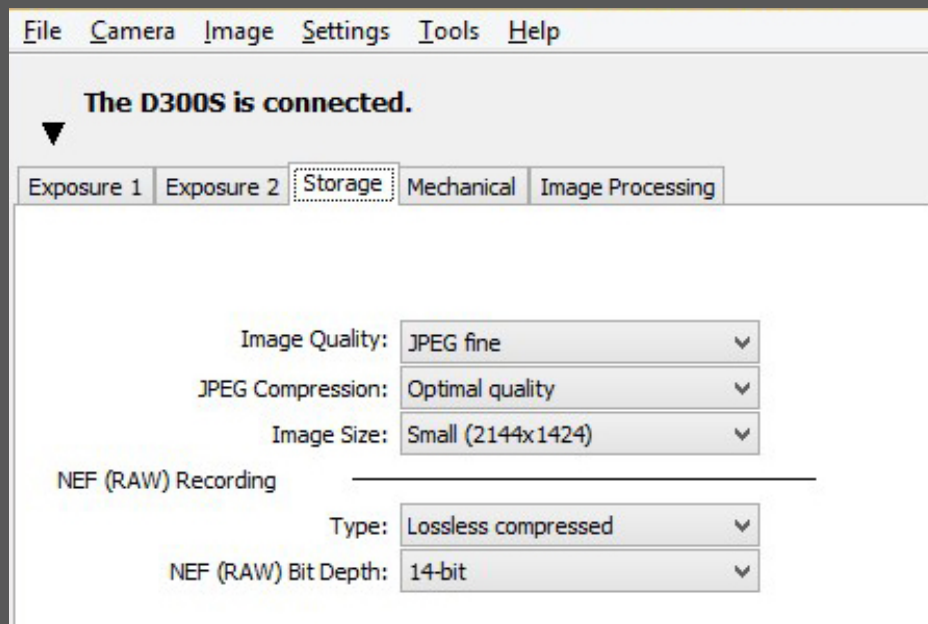

Mechanical:


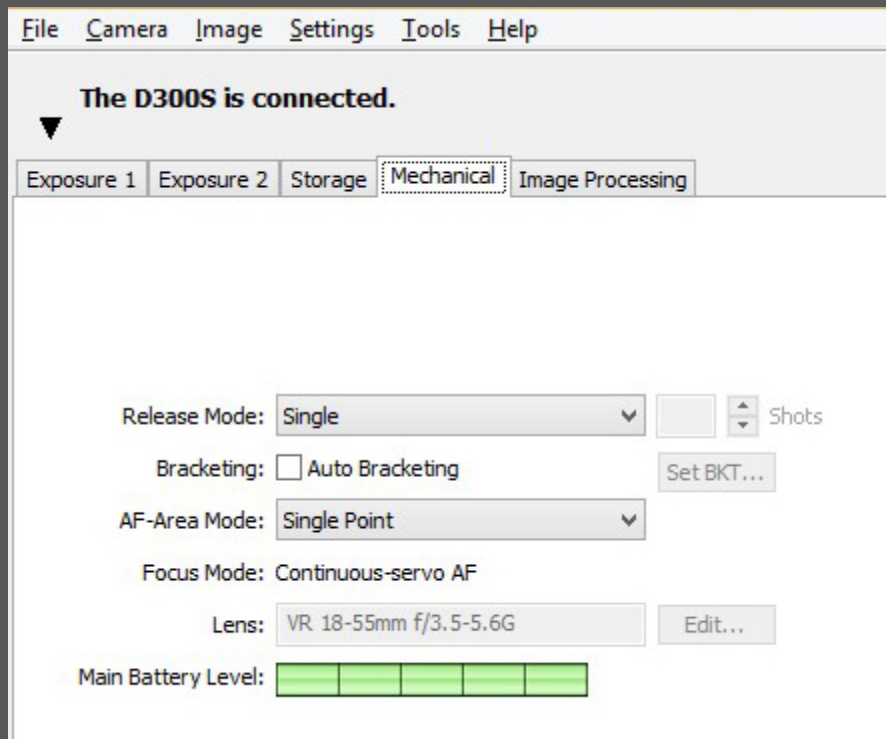


Image Processing:

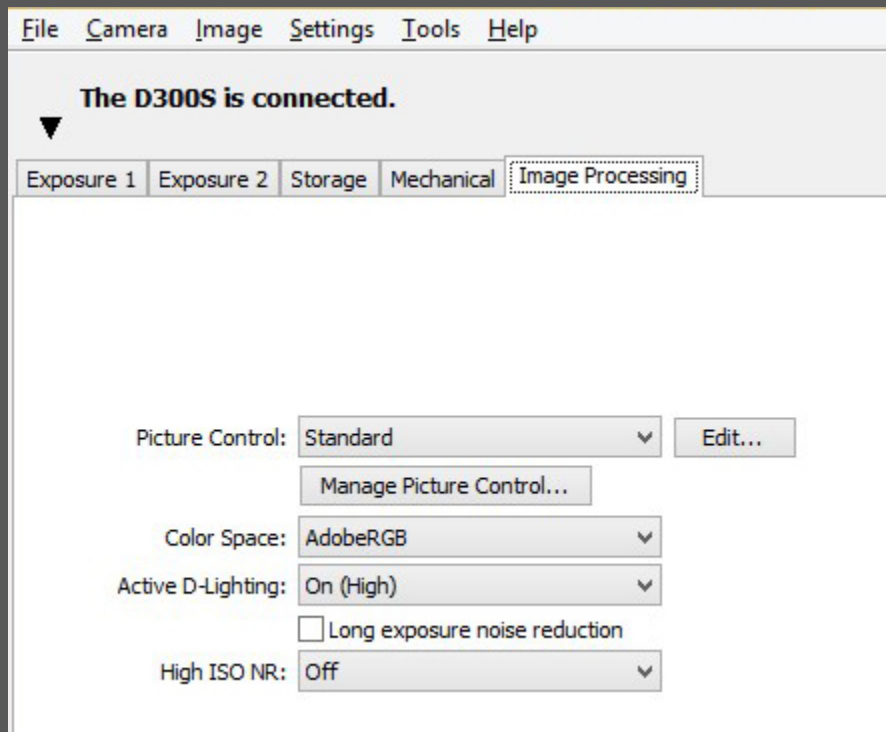


Adjustments for D90


NOTE: It is necessary, that you adjust time and date on your camera manually when bringing it to service for the first time. Otherwise the CCP2 will not detect the camera!

NOTE: Please ensure that you set the photography mode on your camera on „M“ (Manual), so quickPX can take over control over the camera.

1. Please open quickPX at first. (CCP2 will also be started automatically. NOTE: If you are using the trial version of CCP2, please start CCP2 before quickPX!)
2. First please click on the card „Presets“ in the adjustments area on the left of the quickPX interface. You may have to use the arrow-buttons on the right of the card bar to unveil it. Now ensure that the item „Switch live view before shooting“ is checked.

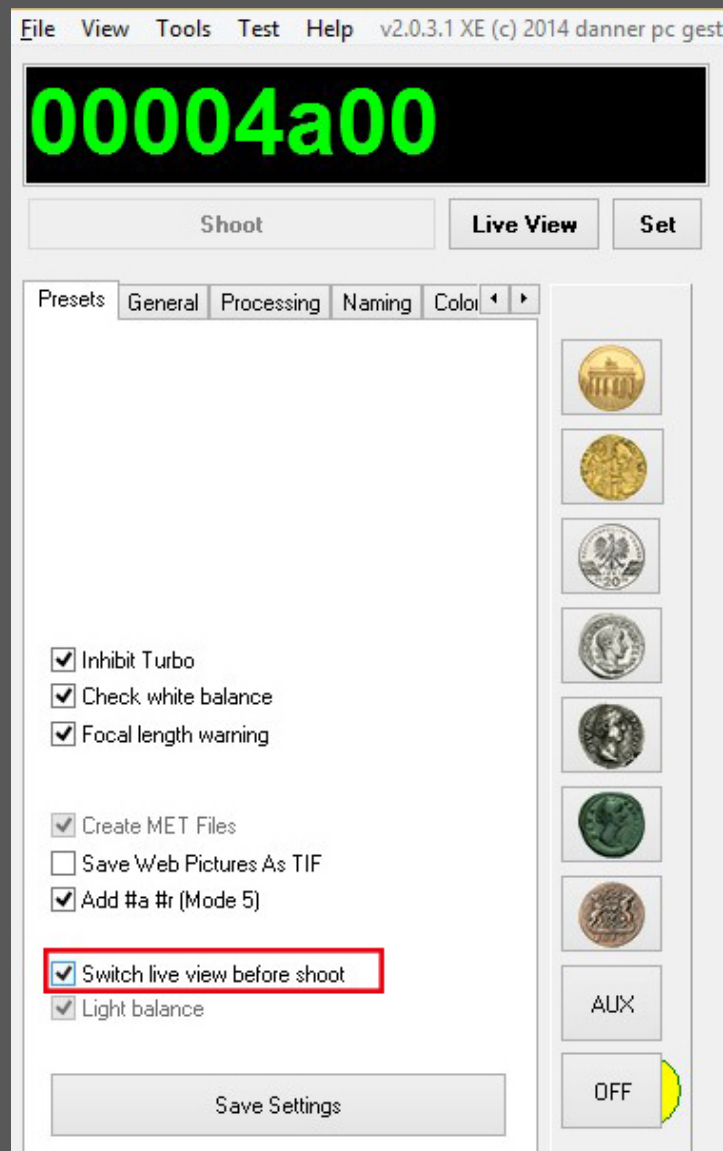


3. To make the remaining adjustments please click on the button „Set“ in the top area of the interface.

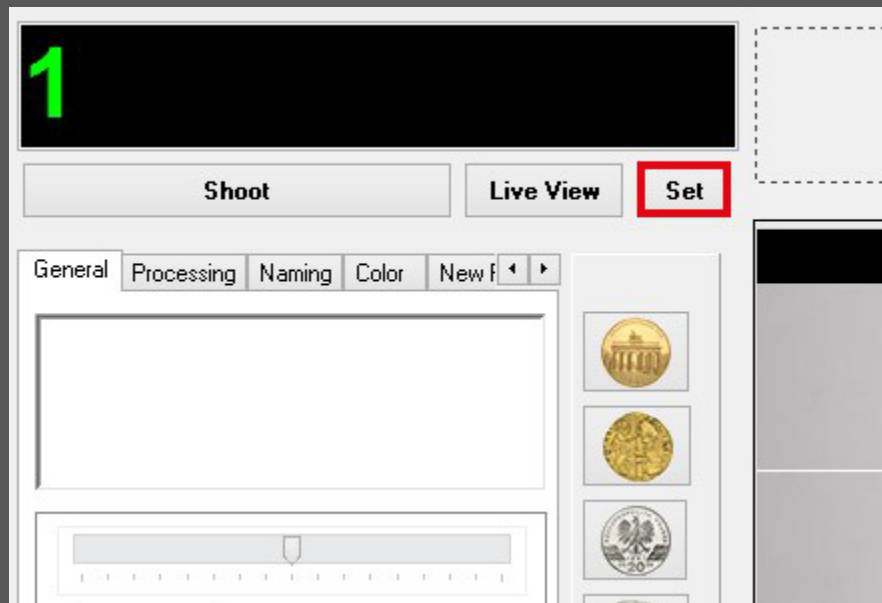


4. Now please ensure that all settings in categories in this menu match the settings depicted in the following pictures.


Exposure 1:

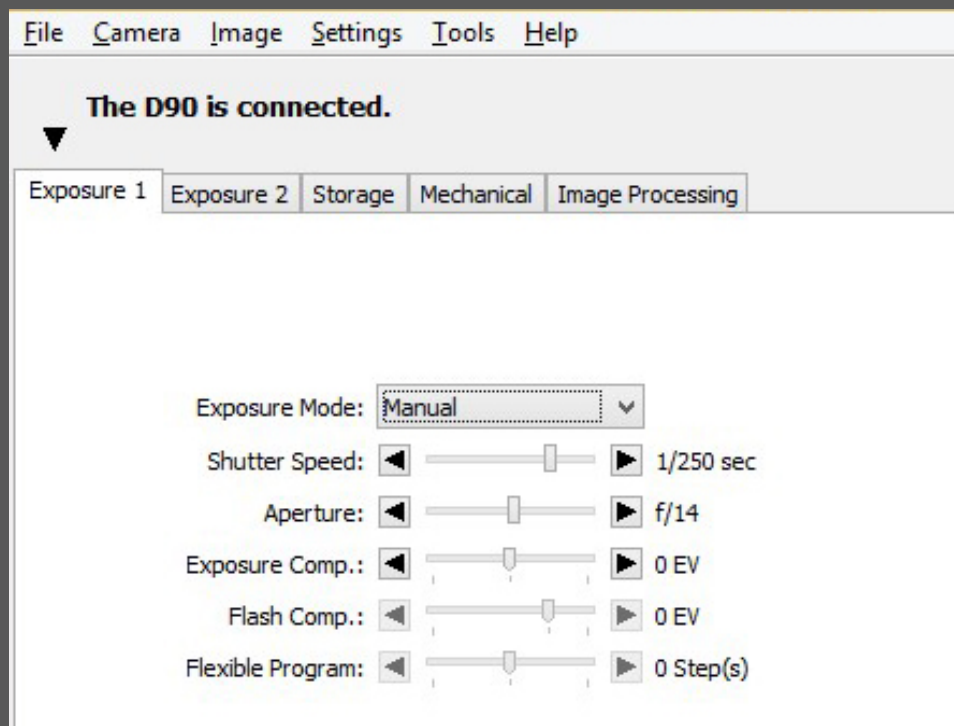

Exposure 2:

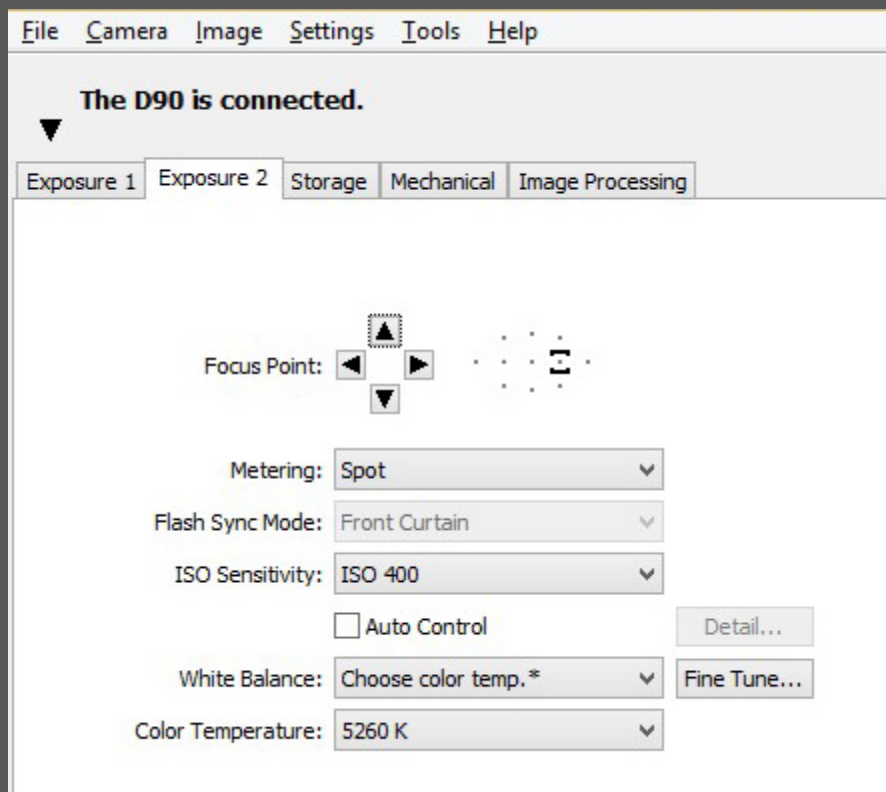

Storage:

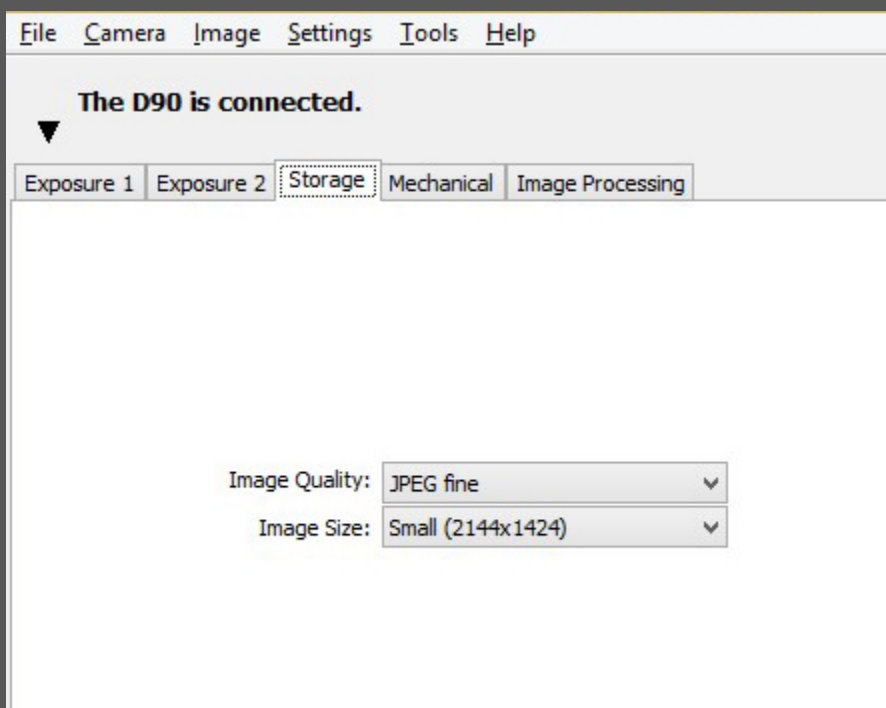

Mechanical:

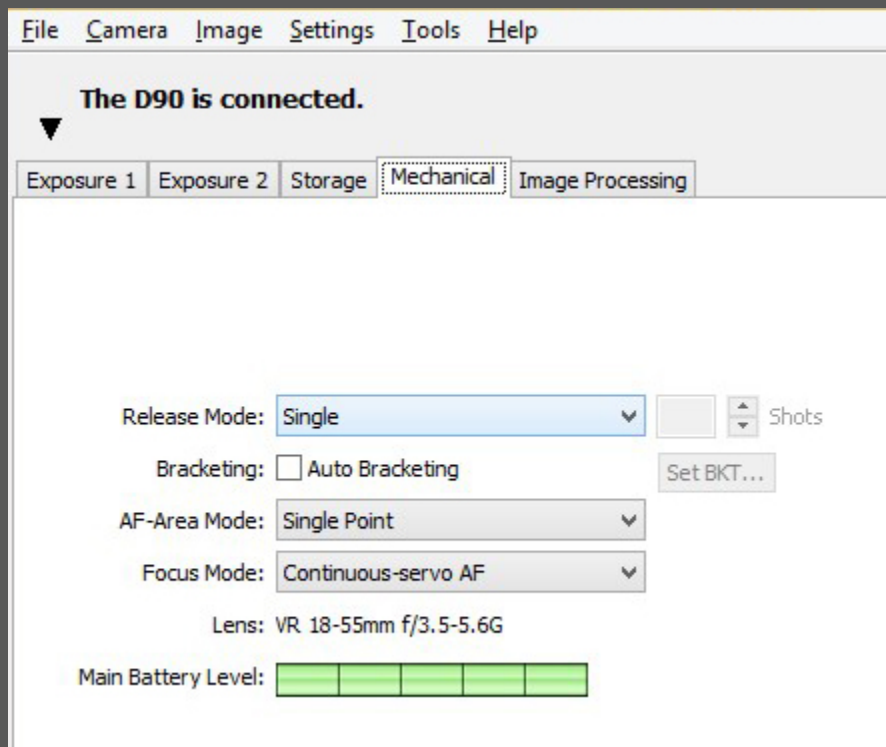


Image Processing:

